

# IP67

## Compact Bench Scale

### INSTRUCTION MANUAL

---

SJ-3000WP/-BT  
SJ-6000WP/-BT  
SJ-15KWP/-BT  
SJ-30KWP/-BT

**AND**

A&D Company, Ltd.

© 2018 A&D Company, Limited. All rights reserved.

No part of this publication may be reproduced, transmitted, transcribed, or translated into any language in any form by any means without the written permission of A&D Company, Limited.

The contents of this manual and the specifications of the instrument covered by this manual are subject to change for improvement without notice.

# CONTENTS

1. SAFETY PRECAUTIONS .....	2
2. PARTS DESCRIPTION .....	3
3. PREPARATION .....	4
3.1. Installing / Exchanging Batteries .....	4
3.2. Setting Up The Scale .....	4
4. DISPLAY AND SYMBOLS.....	5
4.1. Display .....	5
4.2. Symbols .....	5
4.3. Operations And Functions Of Switches .....	6
5. OPERATION .....	7
5.1. Basic Weighing Operation .....	7
5.2. Notes About Operations.....	7
5.3. Weight Display Resolution .....	8
6. SELECTING A WEIGHING UNIT .....	9
6.1. Storing The Weighing Unit .....	9
6.2. Selecting The Weighing Unit.....	9
7. COUNTING MODE.....	10
8. COMPARATOR.....	11
8.1. The Formula To Compare.....	12
8.2. Entering The Comparator Values.....	14
9. AUTO-TARE .....	17
10. WIRELESS COMMUNICATION FUNCTION (SJ-WP-BT model only).....	18
10.1. Wireless Communication Initial Setting.....	18
10.2. SJ-WP-BT Wireless Communication Specification .....	18
10.3. Data output mode (Prt) (refer to “12.FUNCTION SETTINGS “).....	19
10.4. Command Mode .....	20
11. CALIBRATION.....	28
11.1.Calibration Mode.....	28
11.2.Gravity Acceleration Correction .....	30
11.3.Calibration Using A Weight .....	31
11.4.Restoring The Calibration Data To The Factory Set Values .....	32
12. FUNCTION SETTINGS .....	33
12.1.Setting The Parameters .....	33
12.2.Restoring The Function Settings To The Factory Set Values.....	34
12.3.Function List .....	35
13. MAINTENANCE .....	38
13.1.Notes On Maintenance .....	38
13.2. Pick cleaning.....	38
13.3. Error Codes .....	39
14. SPECIFICATIONS.....	40
15. GRAVITY ACCELERATION .....	43

# 1. SAFETY PRECAUTIONS

All safety messages are identified by the following, "WARNING" or "CAUTION", of ANSI Z535.4 (American National Standard Institute: Product Safety Signs and Labels). The meanings are as follows:

 <b>WARNING</b>	A potentially hazardous situation which, if not avoided, could result in death or serious injury.
 <b>CAUTION</b>	A potentially hazardous situation which, if not avoided, may result in minor or moderate injury.

- This manual is subject to change without notice at any time to improve the product.
- Product specifications are subject to change without any obligation on the part of the manufacturer.
- When using the SJ-WP/-BT series, the following safety precautions should always be followed.

## **WARNING**

**Internal service or adjustment to this product should be performed by a qualified person.**

## **CAUTION**

Avoid installing the scale in direct sunlight, which may cause discoloration or malfunctions.

Do not mix battery types, or new and old batteries. Replace with all new batteries at the same time.

If the scale is not to be used for a long period of time, remove all batteries from the battery compartment to avoid leakage.

Avoid overloading the scale.

Avoid using the weighing platform to move the scale, as that could cause damage to the scale.

Avoid chemical solvents. Clean the scale with water.

IP67




Degrees of protection against water: Protected against temporary submergence.

Degrees of protection against solid foreign objects: Dust-tight.

International Protection of IEC60529.

# 2. PARTS DESCRIPTION

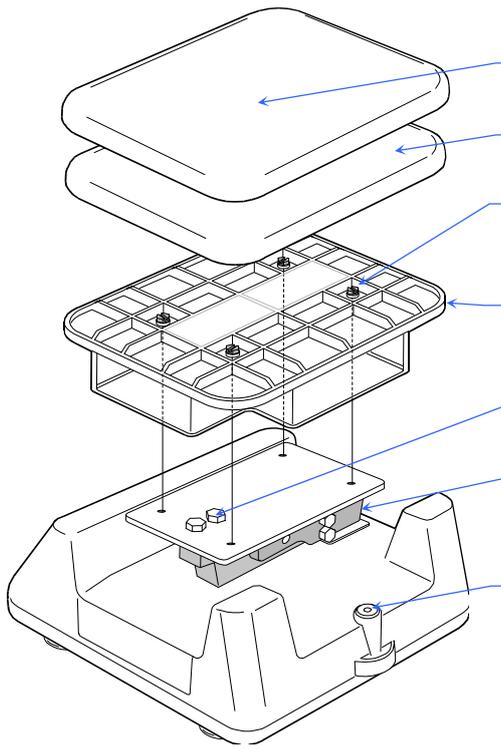

Weighing pan (SUS304)

Inner pan (ABS)

Screw to fix the pan support


Pan support  
The pan support can be removed at times such as when the scale is cleaned.

Do not remove the screws from the scale.

Weight sensor (Load cell)

Bubble Spirit level  
Adjust the leveling feet so that the bubble is in the center of circle.




Switch cover  
When this cover is opened, the **CAL** switches inside are revealed. Keep the cover closed to maintain waterproof performance.

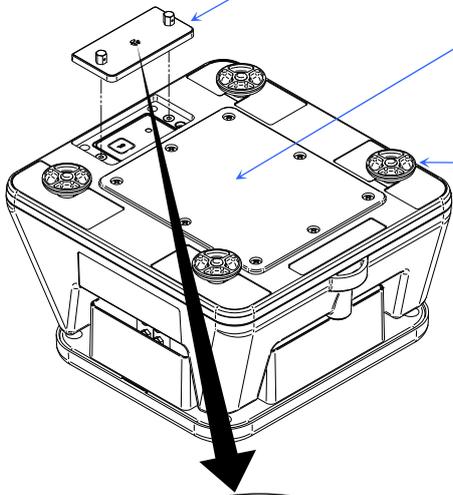

Battery cover  
To maintain waterproof performance, close completely this cover and secure by screws.

Leveling foot


Display and switch

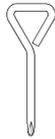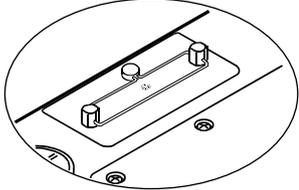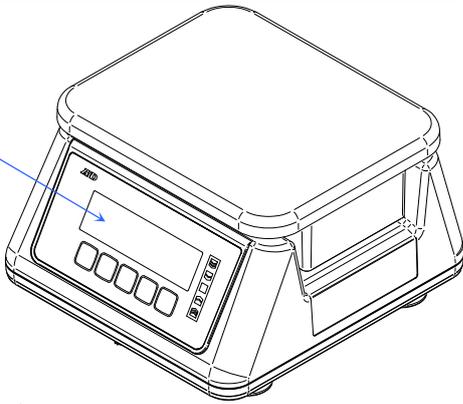

Screwdriver  
Use this screwdriver to open or close the battery cover.

Sealing diagrams

Cleaning picks

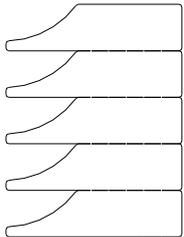


## 3. PREPARATION

### 3.1. Installing / Exchanging Batteries

The batteries are not included with the product. Prepare 6 x "D" size (R20P or LR20) dry-cell batteries.

When the  mark is displayed, exchange the old batteries with new ones.

1. Loosen the screws for securing the battery cover using the provided screwdriver, and remove the battery cover.
2. Remove all the old batteries from the battery compartment.
3. Install six new batteries properly according to the + and - indicators of polarity in the battery compartment.
4. Close the battery cover, and tighten the screws for securing the battery cover.

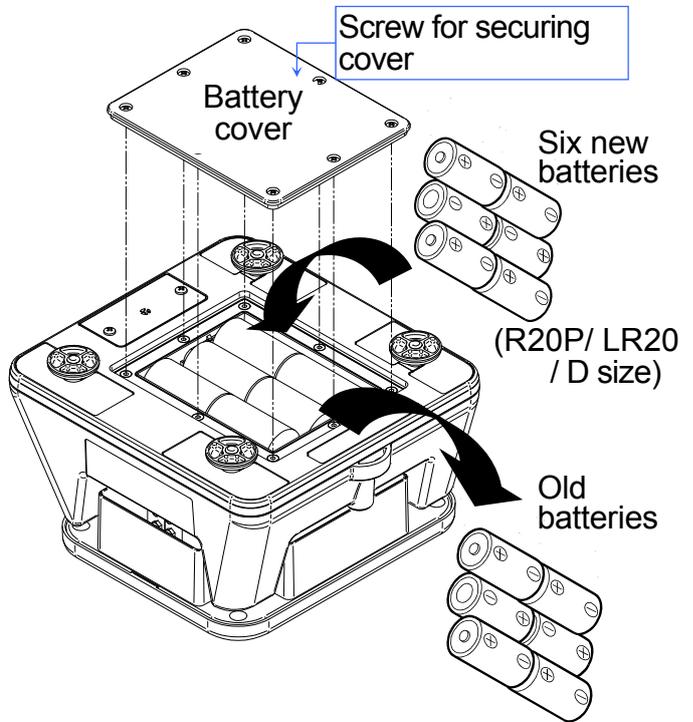

#### CAUTION

- Do not mix used and new batteries. Do not mix the different types of batteries. That may cause damage to the batteries or the scale.
- Take care of the polarity of batteries. The polarity marks are shown in the battery compartment.

### 3.2. Setting Up The Scale

#### CAUTION


- Avoid installing the scale in direct sunlight, that may cause discoloration or malfunctions. Place your SJ-WP/-BT on a firm weighing table so that the scale is level. The scale will not perform accurately when it is not level.
- Place the scale on a firm surface and adjust the feet so that the bubble of the spirit level is in the center of the circle.




## 4. DISPLAY AND SYMBOLS

### 4.1. Display

SJ-WP Model

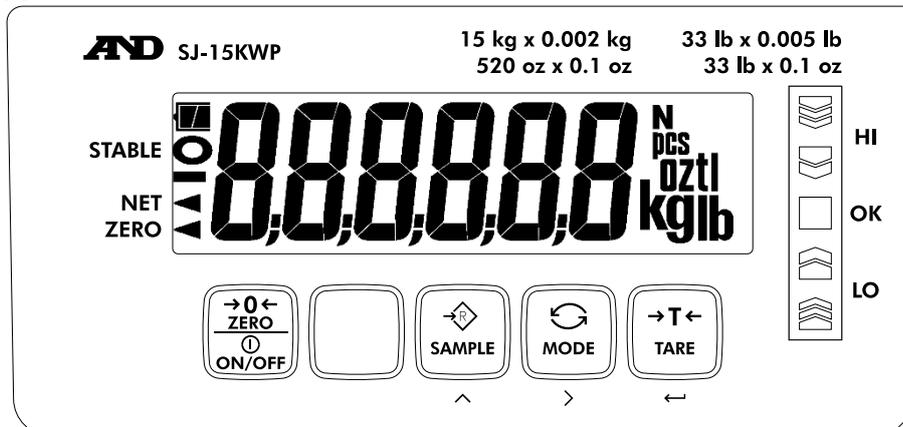



SJ-WP-BT Model

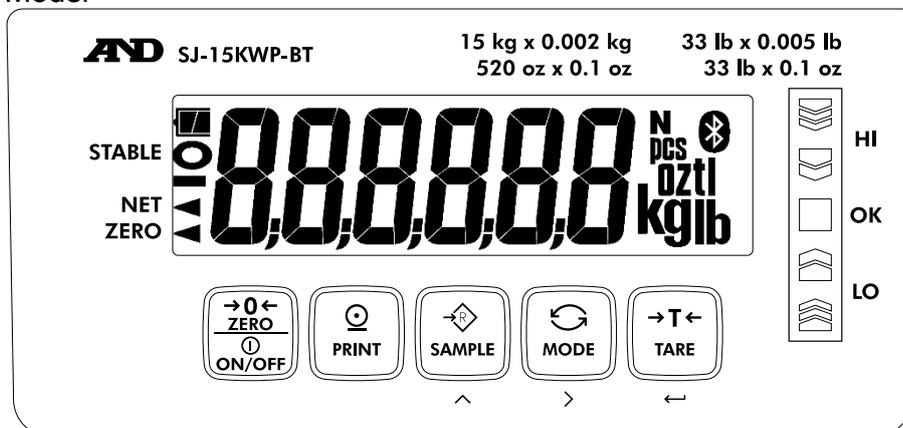

### 4.2. Symbols

Symbol	Description
STABLE	Turns on when the weight value is stable.
NET	Turns on when the NET weight is displayed. (The tare operation is in progress.)
ZERO	Turns on when zero is displayed.
 Comparator LEDs	Turns on when the comparator results are displayed.
Weighing units	"lb", "oz", "ozt", "lb-oz", "tl-s", "tl-h", "tl-t", "t", "pcs", "N", "g" and "kg" are available. A selected unit is displayed.
 Battery indicator	The battery indicator changes as the battery capacity decreases, as shown below:  New → Replace the batteries.
Wireless communication	Turns on when the connection with the wireless communication receiver is established.


### 4.3. Operations And Functions Of Switches

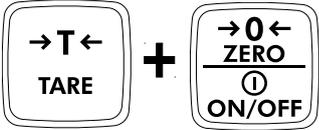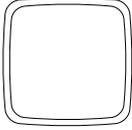Switch	Description
	<p>When the scale is turned <b>off</b>:</p> <p>Press the switch to turn ON the scale. The scale will be automatically set to zero (power-on zero).</p> <p>When the scale is turned <b>on</b>:</p> <p>Press and hold the switch to turn OFF the scale.</p> <p>Press the switch to ZERO the scale and display zero.</p>
	<p>During measurement:</p> <p>In the function setting mode to set the parameters:</p> <p>Press the switch to change the parameter of the selected item.</p>
	<p>During measurement:</p> <p>Press the switch to choose a unit specified in the setting mode.</p> <p>In the function setting mode to set the parameters:</p> <p>Press the switch to select the function item.</p>
	<p>During measurement:</p> <p>Press the switch to tare the scale and display zero (net weight display).</p> <p>In the function setting mode to set the parameters:</p> <p>Other than at the item "unit":</p> <p>Press the switch to store new parameters and return to the weighing mode.</p> <p>At the item "unit":</p> <p>Press the switch to select active / inactive for the displayed unit.</p>
	<p>When the scale is turned <b>off</b>:</p> <p>Press and hold the <span style="border: 1px solid black; padding: 2px;">ON/OFF</span> switch while pressing and holding the <span style="border: 1px solid black; padding: 2px;">TARE</span> switch to enter the function setting mode.</p> <p>Further to the above, continue to press and hold the <span style="border: 1px solid black; padding: 2px;">TARE</span> switch to restore the function settings to the factory set values.</p>
	<p>Output the display value data (SJ-WP-BT model only)</p>
	<p>No use (SJ-WP model only)</p>
	<p>When the scale is turned <b>on</b>:</p> <p>By pressing the switch, the scale proceeds to the calibration mode.</p>

# 5. OPERATION

## 5.1. Basic Weighing Operation

1. Press the **ZERO ON/OFF** switch to turn the power ON.  
All the symbols except the comparator LEDs are displayed. When the weight value becomes stable, the display turns off for a moment and displays zero (power-on zero) with the weighing unit used last before turning off.
2. Select a weighing unit using the **MODE** switch.
- See "6. SELECTING A WEIGHING UNIT" in detail.
3. When the display doesn't show zero, press the **ZERO ON/OFF** switch to set the display to zero.
4. When using a tare (container), place the container on the weighing pan, and press the **TARE** switch to set the display to zero.
5. Place the object to be weighed on the weighing pan or in the container, and wait for the STABLE indicator to turn on and read the value.
6. Remove the object from the weighing pan.
7. Press and hold the **ZERO ON/OFF** switch to turn the power off.

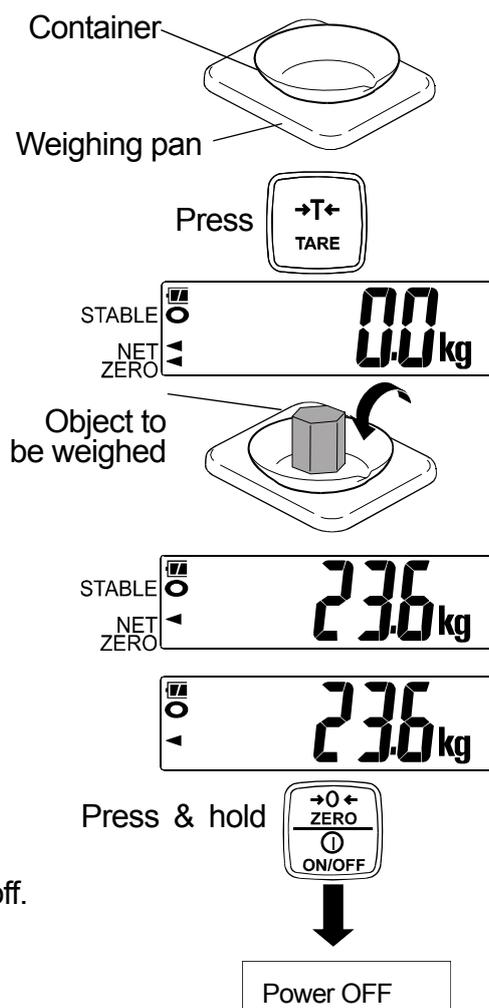

## 5.2. Notes About Operations

### Power-on zero

- If the weight is unstable at power ON, the scale displays **-----**. Check anything touches the weighing pan, or check if there is strong wind or vibration.
- The range for power-on zero is within  $\pm 50\%$  ( $\pm 10\%$  for Legal for Trade models) of the weighing capacity (kg) at the calibrated zero point. If the scale is powered on with a load beyond this range, the scale displays **-----**. Remove the load from the weighing pan.

### ZERO and TARE

- The **ZERO ON/OFF** (as ZERO switch) and **TARE** switches work when the weight value is stable.
- The **ZERO ON/OFF** switch will zero the scale if the weight value is within  $\pm 2\%$  of the weighing capacity (kg) at the power-on zero point. The **ZERO**  $\blacktriangleleft$  indicator turns on. (ZERO operation)

- ❑ The **TARE** switch will tare the scale and subtract the weight to zero as a tare weight when the weight is a plus value. In this case the ZERO ◀ and NET ◀ indicators turn on. (TARE operation) At the zero point, the net weight display shows the tare weight in negative, and the ZERO ◀ and NET ◀ indicators turn on.  
(**Note: In some countries or areas, the ZERO ◀ indicator will not turn on while the scale is tared.**)
- ❑ When the scale is tared, weighing range for net loads is reduced by the amount of the tare weight.
- ❑ When the ZERO operation is performed in the net weight display, the tare operation previously done is cleared and the NET ◀ indicator turns off.  
(**Note: In some countries or areas, the ZERO operation will not clear the TARE operation. Press the **TARE** switch after zeroing the scale with nothing on the weighing pan.**)

### Auto power-off function

- ❑ If no switch is pressed and the stable indicator is displayed for a certain period of time, the scale will automatically turn off. See the function setting **PoFF** to set the elapsed time to turn off.
- ❑ When **E** or **-E** is displayed (refer to "12.3. Error Codes"), the auto power-off function is enabled.

### LCD backlight

- ❑ The LCD backlight is controlled by the functions **L-it** and **L-i**.
- ❑ If no switch is pressed and the weight display continues to be stable for a certain period of time, the LCD backlight will automatically turn off. The elapsed time to turn off is set by the function setting **L-it**. The backlight always on or off is also selectable.
- ❑ The function setting **L-i** adjusts the brightness of the backlight.

## 5.3. Weight Display Resolution

The weight display resolution is a ratio of the minimum display to the weighing capacity. The SJ-WP/-BT series has four types of weight display resolution, as shown below.

Low:	1/3,000
Normal:	1/6,000 or 1/7,500 (depending on the weighing capacity)
High:	1/12,000 or 1/15,000 (depending on the weighing capacity)
Maximum:	1/30,000

The factory setting is the normal resolution. Select the resolution according to your own application in the function setting **rE50**.

- ❑ For details about the minimum display and the weighing capacity, refer to "14. SPECIFICATIONS".
- ❑ The weight display resolution of the Legal for Trade models is fixed. The selection in the function setting **rE50** is not available.
- ❑ In the counting mode, the scale works with the maximum resolution regardless of the weight display resolution selected in the function setting **rE50**.

## 6. SELECTING A WEIGHING UNIT

### 6.1. Storing The Weighing Unit

1. Press the **ZERO ON/OFF** switch while pressing and holding the **TARE** switch in order to display **P-\*\*\*** in the function setting mode.
2. Press twice the **MODE** switch to display **Unit**.
3. Press the **SAMPLE** switch to display a unit. Press the **TARE** switch to activate or deactivate the unit. The indicator **○** is displayed for each active unit.
4. Repeat step 3 for other units.
5. Press the **MODE** switch to move to the next function item when finishing the selection.
6. Press the **TARE** switch to store new units. The scale returns to the weighing mode.

Unit	Symbol	Conversion to gram
Gram	<b>g</b>	1.00000 g
Kilogram	<b>kg</b>	1000.00 g
Pound (UK)	<b>lb</b>	453.59237 g
Ounce (avoir)	<b>oz</b>	28.349523125 g
Troy ounce	<b>ozt</b>	31.1034768 g
Pound - Ounce	<b><sup>oz</sup>lb</b>	—
Tael (Hong Kong general, Singapore)	<b>tl</b> <b>Unit S</b>	37.7994 g
Tael (Hong Kong jewelry)	<b>tl</b> <b>Unit H</b>	37.4290 g
Tael (Taiwan)	<b>tl</b> <b>Unit t</b>	37.5 g
Tola	<b>t</b>	11.6638038 g
Counting unit	<b>pcs</b>	—
Newton	<b>N</b>	See below


Newtons is a value calculated as follows:

$$\text{Newtons} = (\text{value in grams}) \times (9.80665 \text{ m/s}^2) / 1000$$

### 6.2. Selecting The Weighing Unit

In the weighing mode, press the **MODE** switch to select a weighing unit.

Each time the **MODE** switch is pressed, the unit changes as shown below.

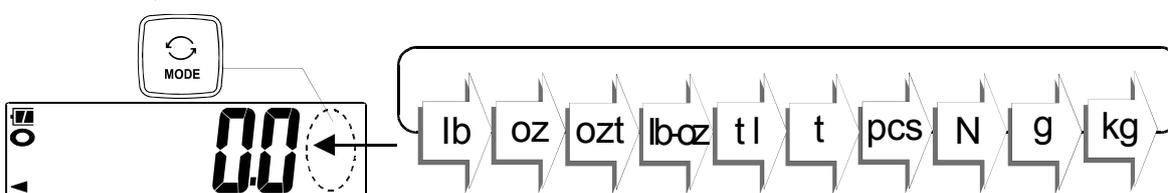


# 7. COUNTING MODE

Determines a unit weight (the weight of one piece) from a known sample quantity, and calculates how many pieces are on the weighing pan using the unit weight. The unit weight is maintained even if the power is turned OFF.

1. Press the **MODE** switch to select "pcs". ("pcs" = pieces)
2. Press the **SAMPLE** switch to enter the sample unit weight storing mode. The numerical value on the left indicates the number of samples.
3. Pressing the **SAMPLE** switch allows you to change the number of samples in the order 5 → 10 → 20 → 50 → 100 → ESC → 5.
  - By pressing the **MODE** switch when **ESC** is displayed, the scale exits the sample unit weight storing mode to the count display.
4. When "-" appears at the right side of the number of samples, press the **ZERO ON/OFF** switch to zero the scale. If necessary, place a container on the weighing pan, and press the **TARE** switch. Confirm that the right side of the number of samples shows zero.
5. Place the correct number of samples on the weighing pan or in the container.
6. Confirm that the STABLE indicator is turned on. Press the **MODE** switch to calculate and store the unit weight. Remove the samples. The scale is set to count objects with this unit weight.
- The total weight of samples should be more than shown below, regardless of the number of samples.
  - SJ-3000WP/-BT: 2.5 g
  - SJ-6000WP/-BT: 5 g
  - SJ-15KWP/-BT: 12.5 g
  - SJ-30KWP/-BT: 25 g

If not, the display shows **Lo ut** and returns to the display of step 5. Increase the number of samples (step 3) and try again.

7. Place the objects to be counted on the weighing pan.

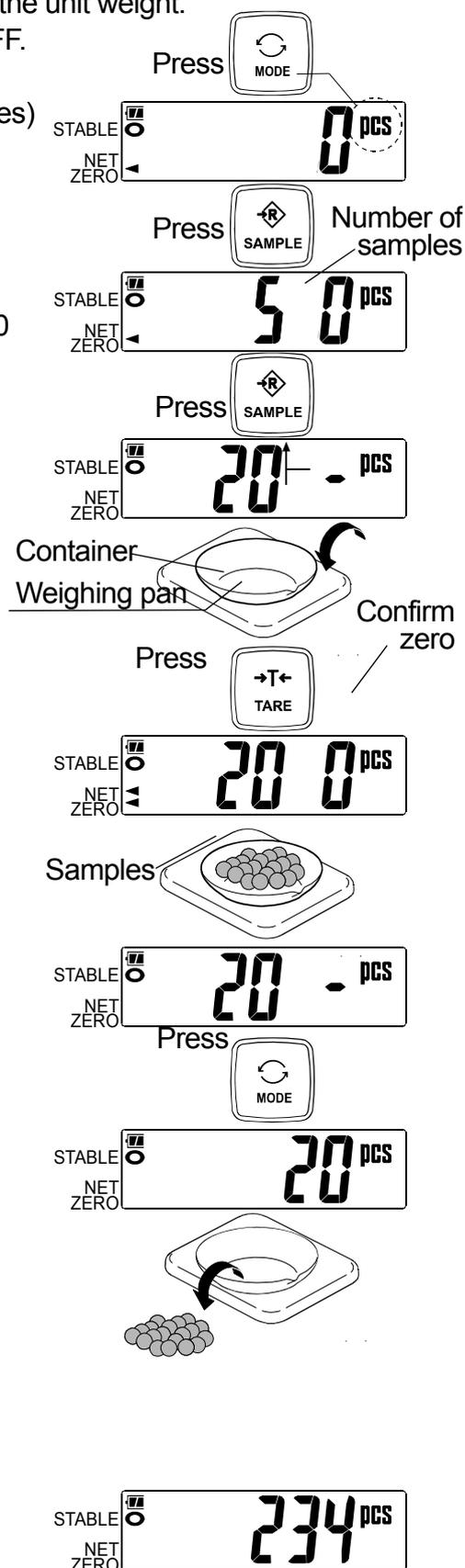


## 8. COMPARATOR

The scale has three-, five- and seven-level comparators.


Each comparator mode compares the weight value against the preset limit values and outputs the results using LEDs (yellow / green / red).

**Note: When the unit is "lb-oz" or "tl", this function can not be used.**

- Five-level comparator mode:  
Uses four comparator values to compare the weight value and outputs results in five levels of LOLO, LO, OK, HI and HIHI.

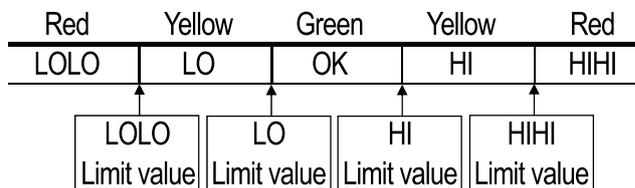


- Three-level comparator mode:  
Uses two comparator values (upper and lower limit values) to compare the weight value and outputs results in three levels of LO, OK and HI.

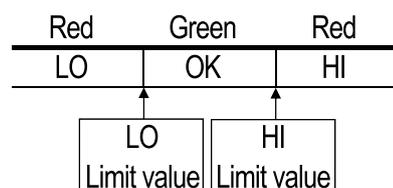

- Seven-level comparator mode (portion weighing mode):

Uses six comparator values to compare the weight value and outputs results in seven levels of outside the lowest limit, level 1 (LOLO), level 2 (LO), Level 3 (OK), level 4 (HI), level 5 (HIHI) and outside the highest limit.

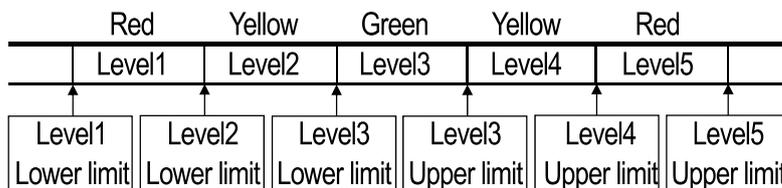

- ❑ To use the comparator modes, the function settings  $\boxed{[P-L]}$  and  $\boxed{[P]}$  must be specified and the comparator values must be set.
- ❑ Using the function setting  $\boxed{[P-L]}$ , select a comparator mode.
  - 0: Five-level comparator mode (Result LED blinks)
  - 1: Five-level comparator mode (Result LED lights)
  - 2: Three-level comparator mode (Result LED blinks)
  - 3: Three-level comparator mode (Result LED lights)
  - 4: Seven-level comparator mode (Result LED blinks)
  - 5: Seven-level comparator mode (Result LED lights)
- ❑ Using the function setting  $\boxed{[P]}$ , select comparison conditions.
  - 0: No comparison (comparator disabled).
  - 1: To compare all data.
  - 2: To compare all stable data.
  - 3: To compare all data which are  $\geq +5d$  or  $\leq -5d$ .
  - 4: To compare stable data which are  $\geq +5d$  or  $\leq -5d$ .
  - 5: To compare all data which are  $\geq +5d$ .
  - 6: To compare stable data which are  $\geq +5d$ .


d = minimum display in kg (Refer to "14. SPECIFICATIONS".)

In the counting mode, "d" is equal to the minimum weight display of kg mode.

## 8.1. The Formula To Compare

Comparison is performed using the formula listed below and the results are output.



### Five-level comparator mode

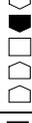Results	Comparison formula	LED display
LOLO	Displayed value < LOLO limit, $\boxed{-\epsilon}$	 ( Red LED on)
LO	LOLO limit $\leq$ Displayed value < LO limit	 ( Yellow LED on)
OK	LO limit $\leq$ Displayed value $\leq$ HI limit	 ( Green LED on)
HI	HI limit < Displayed value $\leq$ HHI limit	 ( Yellow LED on)
HHI	HHI limit < Displayed value, or $\boxed{\epsilon}$	 ( Red LED on)

### Three-level comparator mode

Results	Comparison formula	LED display
LO	Displayed value < LO limit, or $\boxed{-\epsilon}$	 ( Red LED on)
OK	LO limit $\leq$ Displayed value $\leq$ HI limit	 ( Green LED on)
HI	HI limit < Displayed value, or $\boxed{\epsilon}$	 ( Red LED on)

## Seven-level comparator mode (portion weighing mode)

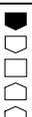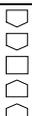Results	Comparison formula	LED display
None	Displayed value < Level 1 lower limit, or $\boxed{-E}$	 ( No LEDs on)
LOLO (Level 1)	Level 1 lower limit $\leq$ Displayed value < Level 2 lower limit	 ( Red LED on)
LO (Level 2)	Level 2 lower limit $\leq$ Displayed value < Level 3 lower limit	 ( Yellow LED on)
OK (Level 3)	Level 3 lower limit $\leq$ Displayed value $\leq$ Level 3 upper limit	 ( Green LED on)
HI (Level 4)	Level 3 upper limit < Displayed value $\leq$ Level 4 upper limit	 ( Yellow LED on)
HIHI (Level 5)	Level 4 upper limit < Displayed value $\leq$ Level 5 upper limit	 ( Red LED on)
None	Level 5 upper limit < Displayed value, or $\boxed{E}$	 ( No LEDs on)

- The comparator values are common to the weighing and counting mode.
- Ignore the decimal point when setting the comparator values.

Example for SJ-6000WP/-BT when the setting value is "001000":

Display mode	Limit value	Capacity / Minimum display
Normal resolution kg	1.000 kg	6.000 kg / 0.001 kg
High resolution kg	0.1000 kg	6.0000 kg / 0.0005 kg
Maximum resolution g	100.0 g	6000.0 g / 0.2 g
Low resolution oz	100.0 oz	210.0 oz / 0.1 oz
Normal resolution oz	10.00 oz	210.00 oz / 0.05 oz
High resolution oz	10.00 oz	210.00 oz / 0.02 oz
Counting mode	1000 pcs	


- The comparator values are maintained even if the power is turned OFF.
- The scale does not judge magnitude relation among the comparator values. Even if the wrong values are set, no error will be shown.

## 8.2. Entering The Comparator Values

### How to Operate

1. Press the **ZERO ON/OFF** switch to put the device in the weighing mode.
2. Press and hold the **MODE** switch to display the currently selected memory number.
3. Each time the **SAMPLE** switch is pressed, the memory number display will be switched.

The currently selected memory number is indicated by the “O” mark being lit.

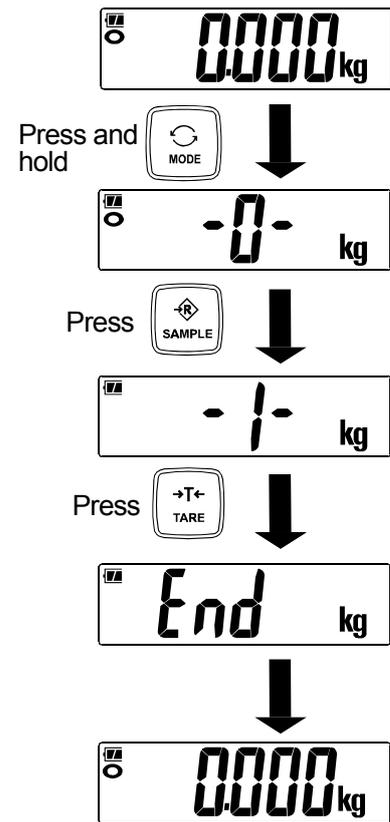

Start comparing using the second memory.


### Selecting the Memory Number

4. Press the **SAMPLE** switch to display the memory number to be changed.
5. Press the **TARE** switch to change the memory number.  
After displaying **End**, the device returns to weighing mode.
6. Start comparing using the memory number changed to.

### Confirming the Upper and Lower Limit Value

4. Press the **SAMPLE** switch to display the memory number to be confirmed.
5. By pressing the **MODE** switch, LO is lit and the lower limit value of the memory number selected is displayed.
6. By pressing the **SAMPLE** switch, HI is lit and the upper limit value of the memory number selected is displayed.
7. To return to the memory number display, press the **SAMPLE** switch.
8. To return to the weighing mode, press the **TARE** switch. (Start comparing using the memory number displayed at this time.)

Confirming the second upper and lower limit value

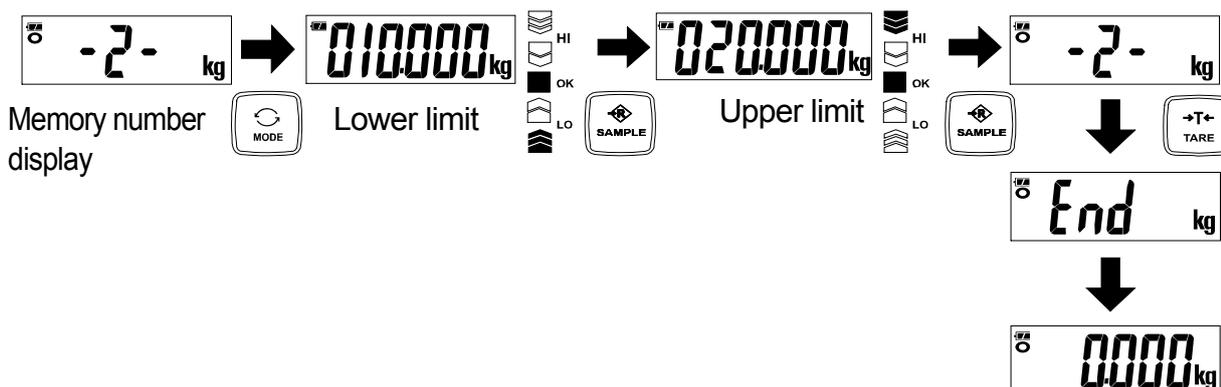

## Setting the Upper and Lower Limit Value

□ When the key lock function active, these operations cannot be used.


4. Press the **SAMPLE** key to display the memory number to be set.
5. Press the **MODE** key to display the lower limit value.
6. Press the **TARE** key at the lower limit value display to make LO and a digit of the value blink.
7. Set the lower limit value by using the following keys.
 

MODE	:	To change which digit is blinking.
SAMPLE	:	To increase by +1 the value of the blinking digit.

The minus sign can be set at the next digit of the least significant digit.  
 The **SAMPLE** switch alternates the minus sign on and off.  
 The blinking "-" shows minus and no sign shows plus.
8. Press the **TARE** key to store the lower limit value. The scale then displays the upper limit value after displaying **End**.
9. Set the upper limit value by using the following keys.
 

MODE	:	To change which digit is blinking.
SAMPLE	:	To increase by +1 the value of the blinking digit.

The minus sign can be set at the next digit of the least significant digit.  
 The **SAMPLE** switch alternates the minus sign on and off.  
 The blinking "-" shows minus and no sign shows plus.
10. Press the **TARE** key to store the upper limit value. The scale then returns to the memory number display.
11. To return to weighing mode, press the **TARE** key. (Start comparing using memory number displayed at this time.)

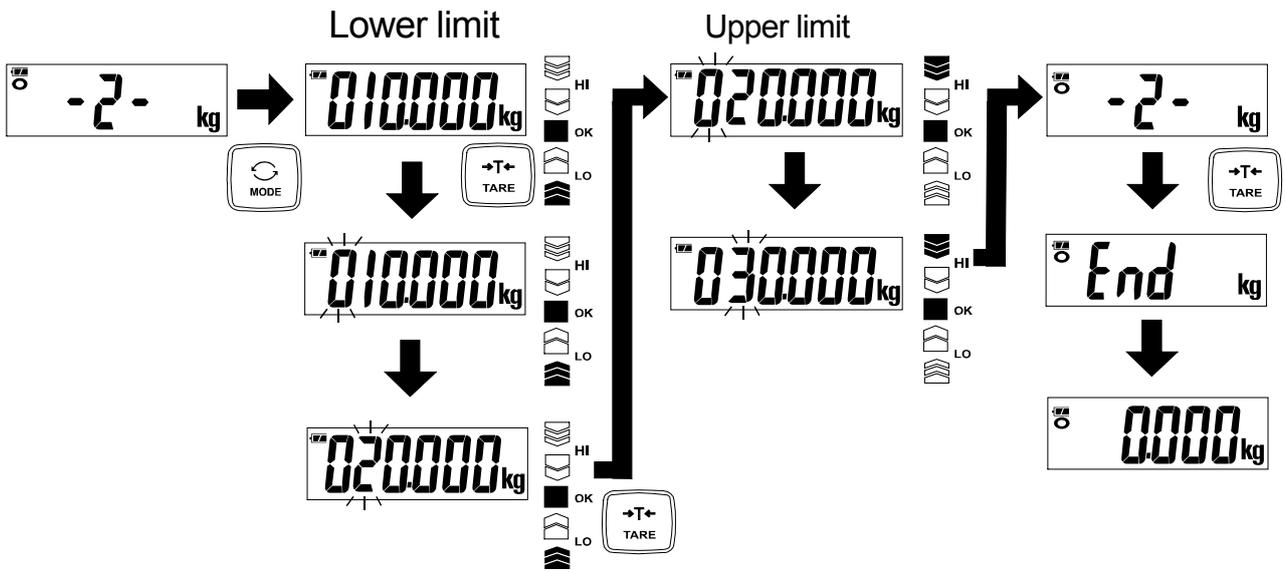

## 9. AUTO-TARE

The SJ-WP/-BT series has an auto-tare function to be used with the comparator mode enabled. If the weight values are in the OK range of comparator limits and stable for a preset period of time, the scale will automatic 1 ally tare the weight and show zero.

- ❑ In some countries or areas, the auto-tare function can not be used on the Legal for Trade models and the selection in the function settings Rt, Rt-t and Rt-F is not available.
- ❑ To use the auto-tare function, set the function settings below.
  - [P 1] : Compare all weighing data (other settings may be used depending on the application).
  - Rt 1 : Auto-tare function enabled.
  - Rt-t 0 to 9 : Select the timing to tare automatically to avoid the wrong tare operation, for example; too early to tare, to take a longer time to go to the next weighing.

### Normal comparison [P-P 0]

- ❑ Start with display zero after tare operation. Place or take away objects until the comparison result will show OK. When the stable indicator is ON for a the preset period of time specified in the function setting Rt-t, the scale will automatically tare the weight, show zero and be ready for next weighing to repeat.

### Negative comparison for take-away [P-P 1] (example with [P-t 0])


- ❑ Take-away check weighing (negative comparison) is the way to compare the negative weight while taking away objects from a container.

Set the function [P-P 1]

together with the auto-tare

function enabled Rt 1. In this operation mode, the scale operates as "take-away the objects" → "OK and stable" → "auto-tare" → "take-away the objects" → ...

In this setting, the polarity of LOLO, LO, HI, and HIHI limit values are ignored and the scale shows the comparator results as below.

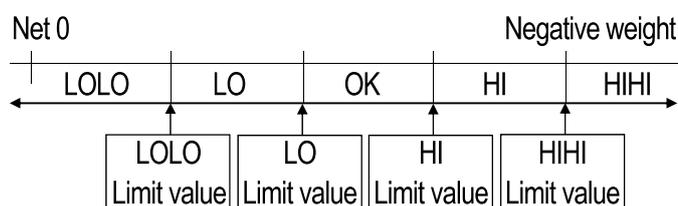

Note: To start the take-away check weighing, be sure to use the TARE switch to tare the weight of the container filled with objects. The ZERO switch may zero the display, and the scale goes below the zero point by taking out the objects. Then, the auto-tare function does not work.

- ❑ When the function "Rt-F 1 Tares the initial (container) weight" is selected: To start the auto-tare function, usually the container (filled with objects) will be placed on the weighing pan and its weight must be tared using the TARE switch. When the function Rt-F 1 is selected, the scale will tare the initial (container) weight automatically. When all load on the weighing pan is removed, the scale will return to the zero point and the tare weight will be automatically cleared. If the scale does not return to the zero point, press the ZERO ON/OFF switch to clear the tare weight.

# 10. WIRELESS COMMUNICATION FUNCTION (SJ-WP-BT model only)

The SJ-WP-BT model has a wireless communication function.

Please purchase optional a wireless communication receiver separately.

This product is not paired with the wireless receiver at the time of shipment. In order to use, it is necessary to perform the procedure in "10.1.Wireless Communication Initial Setting" below.

The SJ-WP-BT model can be paired with one wireless communication device.

## Note

The wireless communication function is a built-in option already incorporated at the time of factory shipment. Therefore, please note that it can not be added to a SJ-WP model later.

## 10.1. Wireless Communication Initial Setting


While the power is off, press and hold the "connection switch" of the wireless communication receiver until the LED lights in orange. (About 3 to 4 seconds)

Turn on the power and wait for a while. If there are multiple scales or balances, turn off the power of those other than the one to be connected. When pairing is successful, the "Wireless communication mark" lights up on the display. If the connection is not successful, please refer to the owner's manual of the wireless communication receiver.

## 10.2. SJ-WP-BT Wireless Communication Specification

**Transmission format** Start-stop synchronous communication(asynchronous), Bidirectional, Half duplex transmission

<b>Signal format</b>	Baud rate	2400bps
	Data bit	7 bit
	Parity	EVEN
	Start bit	1 bit
	Stop bit	1 bit
	Usage code	ASCII
	Terminator	CRLF (CR: 0Dh, LF: 0Ah)

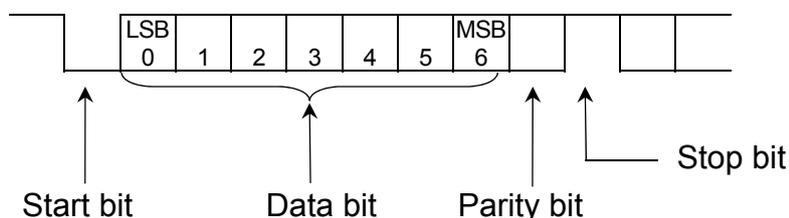

## Data format

S	T	,	+	0	0	0	0	0	.	0	0	_	k	g	C <sub>R</sub>	L <sub>F</sub>
Header		Data										Unit		Terminator		

- There are three kinds of headers for the weighing value as follows.

ST : Weighing data is stable.

US : Weighing data is not stable.

OL : Data is over.(beyond the measurement range)

- Data is always digits including sign, decimal point.

- There are three types of units as follows:

\_ k g : Unit amount of weighing data “kg”

\_ \_ g : Unit amount of weighing data “g”

\_ P C : Unit amount of number “PCS”

- C<sub>R</sub>L<sub>F</sub> is always output for the terminator.

- Example of output data

Weighing data “kg” (+) S T , + 0 0 1 2 . 3 4 5 \_ k g C<sub>R</sub> L<sub>F</sub>

Weighing data “g” (-) S T , - 0 0 0 0 1 2 3 4 \_ \_ g C<sub>R</sub> L<sub>F</sub>

Overweight “kg” (+) O L , + 9 9 9 9 . 9 9 9 \_ k g C<sub>R</sub> L<sub>F</sub>

## 10.3. Data output mode (P<sub>r</sub>t) (refer to “12.FUNCTION SETTINGS”).

### Stream mode (P<sub>r</sub>t - 0)

Data is output continuously. Data output is about 10 times per second.

### Command mode (P<sub>r</sub>t - 0 ~ 5)

The scale is controlled by a command sent from a personal computer etc, connected externally. For details, refer to "10-4. Command mode".

In P<sub>r</sub>t - 1, data is output only by command.

### Output by PRINT key (P<sub>r</sub>t - 2)

When the weighing value is stable (stable mark is lit), pressing the PRINT key will output the data. The display will disappear for a moment to inform you that data has been output.

### Auto print +/- data (P<sub>r</sub>t - 3)

Data is output when the weighing value is stabilized (stability mark is lit) and its value is +5d or greater, or -5d or smaller. The next output will be after the weighing value has returned to the range of -4d to +4d.

#### Auto print + data output ( $P_{rL} - 4$ )

When the value is stabilized (stability mark is lit) and its value is + 5d (d=scale) or greater, data will be output. The next output will be after the weighing value returns to +4d or less.

#### Auto print +/- data and comparison result OK ( $P_{rL} - 5$ )

Data is output when the weighing value is stabilized (stability mark is lit), its value is +5d or greater or -5d or smaller and the comparison result is OK.

The next output will be after the weighing value has returned to the range of -4d to +4d.

#### Auto print + data and comparison result OK ( $P_{rL} - 6$ )

Data is output when the weighing value is stabilized (stability mark is lit), its value is +5d (d=weight minimum display) or greater and the comparison result is OK. The next output will be after the weighing value returns to +4d or smaller.

## 10.4. Command Mode

In the command mode, the scale is controlled by commands that come from an external device such as a computer.

# Command List

Command	Description	Remarks
Q	Requests data be output immediately.	
Z	Zeros the scale when the weighing value is stable.	Same as the <input type="button" value="ZERO"/> key.
T	Tares the scale when the weighing value is stable.	Same as the <input type="button" value="TARE"/> key.
U	Switches the weighing unit.	Same as the <input type="button" value="MODE"/> key.
CT	Clears tare	
?H3	In five-level comparator mode : Not used	The output of setting values for comparator mode *Internal setting in comparator comparison mode  Five-level mode: "[P-L 0]" Three-level mode: "[P-L 2]" Seven-level mode: "[P-L 4]" "[P-L 5]"
	In three-level comparator mode : Not used	
	In seven-level comparator mode : Threshold value of rank 5 is output.	
?H2	In five-level comparator mode : HIHI limit value is output.	
	In three-level comparator mode : Not used	
	In seven-level comparator mode : Threshold value of rank 4 is output.	
?H1	In five-level comparator mode : HI limit value is output.	
	In three-level comparator mode : HI limit value is output.	
	In seven-level comparator mode : Upper threshold value of rank 3 is output.	
?L1	In five-level comparator mode : LO limit value is output.	
	In three-level comparator mode : LO limit value is output.	
	In seven-level comparator mode : Threshold value of rank 3 is output.	
?L2	In five-level comparator mode : LOLO limit value is output.	
	In three-level comparator mode : Not used	
	In seven-level comparator mode : Threshold value of rank 2 is output.	
?L3	In five-level comparator mode : Not used	
	In three-level comparator mode : Not used	
	In seven-level comparator mode : Threshold value of rank 1 is output.	

Command	Description	Remarks
H3	In five-level comparator mode : Not used	Input the six-digit value excluding the polarity and decimal point.
	In three-level comparator mode : Not used	
	In seven-level comparator mode : The threshold value of rank 5 is stored.	
H2	In five-level comparator mode : HHI limit value is stored.	
	In three-level comparator mode : Not used	
	In seven-level comparator mode :The threshold value of rank 4 is stored.	
H1	In five-level comparator mode : HI limit value is stored.	
	In three-level comparator mode : HI limit value is stored.	
	In seven-level comparator mode : The upper threshold value of rank 3 is stored.	
L1	In five-level comparator mode : LO limit value is stored.	
	In three-level comparator mode : LO limit value is stored.	
	In seven-level comparator mode : The lower threshold value of rank 3 is stored.	
L2	In five-level comparator mode : LOLO limit value is stored.	
	In three-level comparator mode :Not used	
	In seven-level comparator mode : The threshold value of rank 2 is stored.	
L3	In five-level comparator mode : Not used	
	In three-level comparator mode : Not used	
	In seven-level comparator mode : The threshold value of rank 1 is stored.	

## Command Examples

("\_" stands for "space"(20H))

### □ To request weighing data

Command 

Q	CR	LF
---	----	----

Reply 

S	T	,	+	0	0	1	2	.	3	4	5	_	k	g	CR	LF
U	S	,	+	0	0	0	7	.	8	9	0	_	k	g	CR	LF
O	L	,	+	9	9	9	9	.	9	9	9	_	k	g	CR	LF

 Stable positive data  
 Unstable positive data  

E
---

 display

### □ To set zero point

Command 

Z	CR	LF
---	----	----

Reply 

Z	CR	LF
---	----	----

 When zero operation is possible

### □ To tare the weighing value

Command 

T	CR	LF
---	----	----

Reply 

T	CR	LF
---	----	----

 When the tare operation can be performed

### □ To cancel tare value

Command 

C	T	CR	LF
---	---	----	----

Reply 

C	T	CR	LF
---	---	----	----

 Clear tare value (including when there is no tare)

### □ In 5-level comparator mode...Not use

In 3-level comparator mode...Not use

In 7-level comparator mode...Outputs upper threshold value of rank 5 in use

Command 

?	H	3	CR	LF
---	---	---	----	----

Reply 

H	3	,	+	0	0	5	0	0	0	CR	LF
---	---	---	---	---	---	---	---	---	---	----	----

### □ In 5-level comparator mode...Outputs the HHI threshold value (upper limit value) in use

In 3-level comparator mode...Not use

In 7-level comparator mode...Outputs upper threshold value of rank 4 in use

Command 

?	H	2	CR	LF
---	---	---	----	----

Reply 

H	2	,	+	0	0	4	0	0	0	CR	LF
---	---	---	---	---	---	---	---	---	---	----	----

- ❑ In 5-level comparator mode...Outputs the HI threshold value (upper limit value) in use
- In 3-level comparator mode...Outputs the HI threshold value (upper limit value) in use
- In 7-level comparator mode...Outputs upper threshold value of rank 3 in use

Command 

?	H	1	C <sub>R</sub>	L <sub>F</sub>
---	---	---	----------------	----------------

Reply 

H	1	,	+	0	0	3	0	0	0	C <sub>R</sub>	L <sub>F</sub>
---	---	---	---	---	---	---	---	---	---	----------------	----------------

- ❑ In 5-level comparator mode...Outputs the LO threshold value (lower limit value) in use
- In 3-level comparator mode...Outputs the LO threshold value (lower limit value) in use
- In 7-level comparator mode...Outputs lower threshold value of rank 3 in use

Command 

?	L	1	C <sub>R</sub>	L <sub>F</sub>
---	---	---	----------------	----------------

Reply 

L	1	,	+	0	0	2	0	0	0	C <sub>R</sub>	L <sub>F</sub>
---	---	---	---	---	---	---	---	---	---	----------------	----------------

- ❑ In 5-level comparator mode...Outputs the LOLO threshold value (lower limit value) in use
- In 3-level comparator mode...Not use
- In 7-level comparator mode...Outputs lower threshold value of rank 2 in use

Command 

?	L	2	C <sub>R</sub>	L <sub>F</sub>
---	---	---	----------------	----------------

Reply 

L	2	,	+	0	0	1	0	0	0	C <sub>R</sub>	L <sub>F</sub>
---	---	---	---	---	---	---	---	---	---	----------------	----------------

- ❑ In 5-level comparator mode...Not use
- In 3-level comparator mode...Not use
- In 7-level comparator mode...Outputs lower threshold value of rank 1 in use

Command 

?	L	3	C <sub>R</sub>	L <sub>F</sub>
---	---	---	----------------	----------------

Reply 

L	3	,	+	0	0	0	0	0	0	C <sub>R</sub>	L <sub>F</sub>
---	---	---	---	---	---	---	---	---	---	----------------	----------------

- ❑ In 5-level comparator mode...Not use
- In 3-level comparator mode...Not use
- In 7-level comparator mode...Outputs upper threshold value of rank 3 in use

Command 

H	3	,	+	0	0	5	0	0	0	C <sub>R</sub>	L <sub>F</sub>
---	---	---	---	---	---	---	---	---	---	----------------	----------------

Reply 

H	3	,	+	0	0	5	0	0	0	C <sub>R</sub>	L <sub>F</sub>
---	---	---	---	---	---	---	---	---	---	----------------	----------------

- ❑ In 5-level comparator mode... Sets the HHI threshold value (upper limit value) in use

In 3-level comparator mode... Not use

In 7-level comparator mode... Sets upper threshold value of rank 4 in use

Command 

H	2	,	+	0	0	4	0	0	0	C <sub>R</sub>	L <sub>F</sub>
---	---	---	---	---	---	---	---	---	---	----------------	----------------

Reply 

H	2	,	+	0	0	4	0	0	0	C <sub>R</sub>	L <sub>F</sub>
---	---	---	---	---	---	---	---	---	---	----------------	----------------

- ❑ In 5-level comparator mode... Sets the HI threshold value (upper limit value) in use

In 3-level comparator mode... Sets the HI threshold value (upper limit value) in use

In 7-level comparator mode... Sets upper threshold value of rank 3 in use

Command 

H	1	,	+	0	0	3	0	0	0	C <sub>R</sub>	L <sub>F</sub>
---	---	---	---	---	---	---	---	---	---	----------------	----------------

Reply 

H	1	,	+	0	0	3	0	0	0	C <sub>R</sub>	L <sub>F</sub>
---	---	---	---	---	---	---	---	---	---	----------------	----------------

- ❑ In 5-level comparator mode... Sets the LO threshold value (lower limit value) in use

In 3-level comparator mode... Sets the LO threshold value (lower limit value) in use

In 7-level comparator mode... Sets lower threshold value of rank 3 in use

Command 

L	1	,	+	0	0	2	0	0	0	C <sub>R</sub>	L <sub>F</sub>
---	---	---	---	---	---	---	---	---	---	----------------	----------------

Reply 

L	1	,	+	0	0	2	0	0	0	C <sub>R</sub>	L <sub>F</sub>
---	---	---	---	---	---	---	---	---	---	----------------	----------------

- ❑ In 5-level comparator mode... Sets the LOLO threshold value (lower limit value) in use

In 3-level comparator mode... Not use

In 7-level comparator mode... Sets lower threshold value of rank 2 in use

Command 

L	2	,	+	0	0	1	0	0	0	C <sub>R</sub>	L <sub>F</sub>
---	---	---	---	---	---	---	---	---	---	----------------	----------------

Reply 

L	2	,	+	0	0	1	0	0	0	C <sub>R</sub>	L <sub>F</sub>
---	---	---	---	---	---	---	---	---	---	----------------	----------------

- ❑ In 5-level comparator mode...Not use
- In 3-level comparator mode...Not use
- In 7-level comparator mode...Sets lower threshold value of rank 1 in use

Command	L	3	,	+	0	0	0	0	0	0	0	0	C <sub>R</sub>	L <sub>F</sub>
Reply	L	3	,	+	0	0	0	0	0	0	0	0	C <sub>R</sub>	L <sub>F</sub>

# Precautions related to radio waves

## FCC

---

Contains Transmitter Module FCC ID: RYYEYSHCN

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation

### FCC WARNING

Changes or modification not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference, in which case the user will be required to correct the interference at his own expense.

## IC

---

### IC RADIATION EXPOSURE STATEMENT FOR CANADA

This device complies with Industry Canada license-exempt RSS standards. Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept any interference, including interference that may cause undesired operation of the device.

L'utilisation de ce dispositif est autorisée seulement aux conditions suivantes:

(1) il ne doit pas produire de brouillage et (2) l'utilisateur du dispositif doit être prêt à accepter tout brouillage radioélectrique reçu, même si ce brouillage est susceptible de compromettre le fonctionnement du dispositif.

This product is certified as type of portable device within Industry Canada Rules. To maintain compliance with RF Exposure requirements, please use within the specification of this product. Ce produit est certifié comme type de l'appareil portable avec Industrie Règles de Canada. Pour maintenir l'acquiescement avec exigence Exposition de RF, veuillez utiliser dans spécification de ce produit.

Contains Transmitter module IC : 4389B-EYSHCN

# 11. CALIBRATION

Adjusts the scale for accurate weighing. Calibrate the scale in the following cases.

- When the scale is first installed.
- When the scale has been moved.
- When the ambient environment has changed.
- For regular calibration.

**Note: The Legal for Trade models can not be re-calibrated if they have been sealed.**

## 11.1. Calibration Mode

- The calibration mode has the following three functions.
  - Gravity acceleration correction
  - Calibration using a weight
  - Restoring the factory set values
- How to enter the calibration mode


### Method 1:


1. Make sure that the scale is in the weighing mode.
2. Press and hold the  SAMPLE  switch until the "CAL" appears and release the switch.

**Note: The above operation is disabled for the Legal for Trade models.**

### Method 2:

1. Make sure that the scale is in the weighing mode.
2. Loosen the two screws on the switch cover and open the switch cover. The calibration (CAL) switch is located inside.
3. Press the  CAL  switch. The scale displays the "CAL".

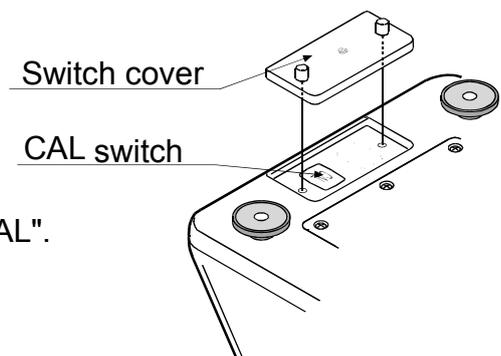

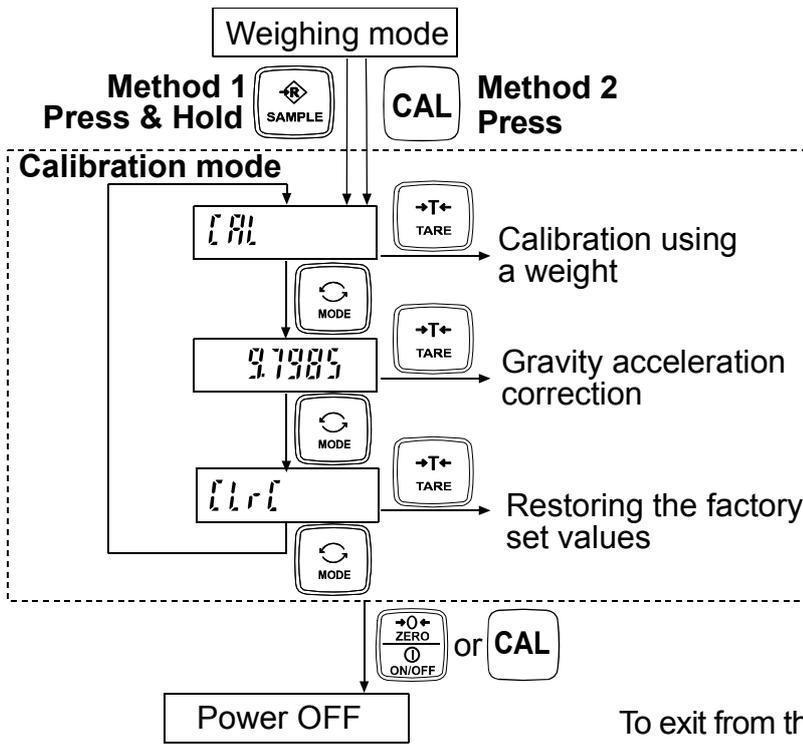

Press the **MODE** switch to select an item

and

Press the **TARE** switch to execute it.


To exit from the calibration mode, press and hold the **ZERO ON/OFF** switch or press the **CAL** switch. The scale turns the power OFF.

## 11.2.Gravity Acceleration Correction

When the scale is first used or has been moved to another location, it should be calibrated using a calibration weight. But if a calibration weight is not available, the gravity acceleration correction will compensate the scale. Change the gravity acceleration value stored in the scale to the value of the area where the scale will be used. Refer to the gravity acceleration map at the end of this manual.

**Note: Gravity acceleration correction is not required when the scale will be calibrated using a calibration weight at the place where it is to be used.**


1. Refer to "Calibration Mode" to enter the calibration mode. The **CAL** is displayed.
2. Press the **TARE** switch to enter the gravity acceleration value setting mode.
3. Change the displayed value using the following switches.  
**MODE** To shift the blinking digit to the right.  
**SAMPLE** To increase the value of the blinking digit by one.
4. Press the **TARE** switch. The display shows **End** and returns to the newly stored gravity acceleration value.
5. When calibration using a calibration weight is to be performed, go to step 3 of "11.3.Calibration Using A Weight".  
To finish the setting procedure, press and hold the **ZERO ON/OFF** switch or the **CAL** switch. The scale returns to the weighing mode.

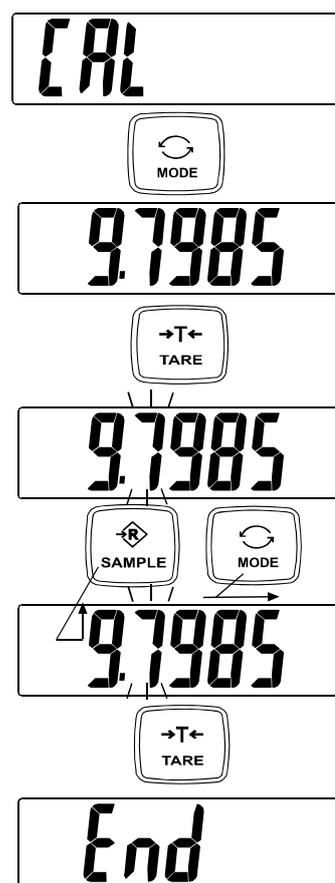

## 11.3. Calibration Using A Weight

Prepare a weight, preferably a weight with the same value as the weighing capacity of the scale to be calibrated. Note that the calibration weight value can be changed.

1. Turn the power ON and warm up the scale for at least half an hour.
- ❑ Change the function setting **[PoFF]** or place something on the weighing pan to disable the auto power-off function.
2. Refer to "Calibration Mode" to enter the calibration mode.  
The **[CAL]** is displayed.
3. Press the **[TARE]** switch, then **[°CAL 0]** is displayed.  
Confirm that nothing is placed on the weighing pan and wait for the STABLE indicator to turn on.
4. Press the **[TARE]** switch. The scale calibrates the zero point and displays the value of the calibration weight (SPAN calibration).
  - ❑ The calibration weight value is equal to the weighing capacity. (factory setting)
  - ❑ If SPAN calibration is not to be performed, turn the power OFF to exit from the calibration procedure.
5. To calibrate with a weight different from the weighing capacity, change the displayed value using the following switches.
  - [MODE]** To shift the digit that is blinking to the right.
  - [SAMPLE]** To increase the value of the blinking digit by one.
- ❑ Using a weight with the same value as the weighing capacity is recommended. If other weights are used, use one with a value greater than two-thirds of the capacity.
6. Place the calibration weight with the same value as displayed on the weighing pan, and wait for the STABLE indicator to turn on.
7. Press the **[TARE]** switch. The scale calibrates SPAN and **[End]** is displayed. Then, the display returns to **[CAL]**.  
To finish the procedure, press and hold the **[ZERO ON/OFF]** switch or press the **[CAL]** switch. The scale turns the power OFF.

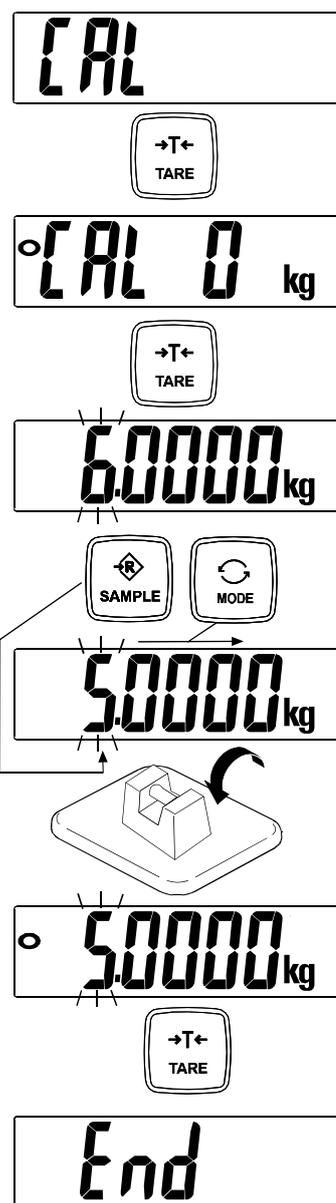

**Note:** If the scale will be moved to another location, set the gravity acceleration value for the present location first and calibrate the scale using a weight. Then, change the gravity acceleration value for the new location.

## 11.4. Restoring The Calibration Data To The Factory Set Values

If the gravity acceleration value or calibration data is changed unintentionally, restore those values to the factory set values, as follows.

1. Refer to "Calibration Mode" to enter the calibration mode. The **CAL** is displayed.
2. Press the **MODE** switch twice to display **LLrL**.
3. Press the **TARE** switch to display **LLrLno** with "no" blinking.
4. Press the **SAMPLE** switch. **LLrLno** changes to **LLrLGo** with "Go" blinking.
- To cancel the restoring procedure, press the **ZERO ON/OFF** switch. The display returns to step 2.
5. When **LLrLGo** is displayed, press the **TARE** switch. The factory set values are restored and **End** is displayed. Then, the display returns to **LLrL**.

To finish the setting procedure, press and hold the **ZERO ON/OFF** switch or press the **CAL** switch. The scale turns the power OFF.

**CAL**

Press twice 

**LLrL**

Press 

**LLrLno**

Press 

**LLrLGo**

Press 


**End**

↓  
**LLrL**

# 12. FUNCTION SETTINGS


The scale has function settings to specify the scale performance.

The parameters set in the function settings are maintained even if the power is turned OFF.

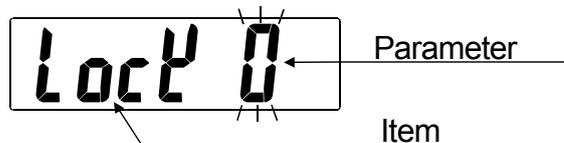


## 12.1. Setting The Parameters

1. Turn the power OFF.
2. **Press** and hold the **TARE** switch and **press** the **ZERO ON/OFF** switch to turn the power ON.  
The software version is displayed.
  - ❑ "\*\*\*" indicates the software version number.
3. After about one second, the item is displayed.
4. Change the item or parameter using the following switches.
  - MODE** To display the next item.
  - SAMPLE** To increase the value of the blinking digit by one (to change the parameter).
- ❑ When the parameter is changed, the STABLE indicator turns off.
5. **Press** the **TARE** switch to store the setting value.  
After displaying **End**, the scale goes to the weighing mode.
  - ❑ To cancel the setting procedure without storing the value, **press** and hold the **ZERO ON/OFF** switch to turn off the scale.

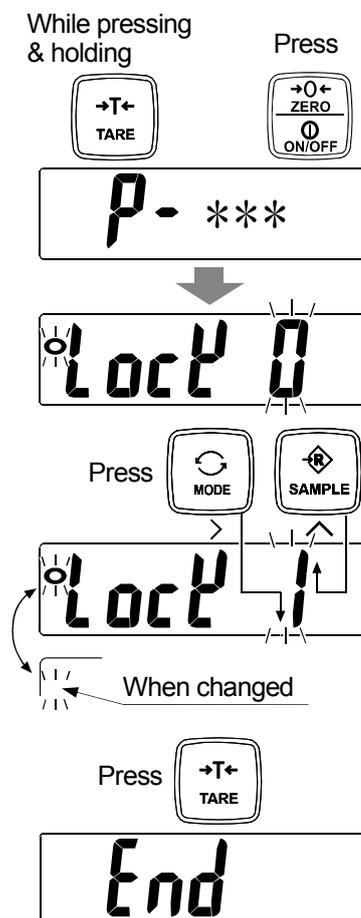

## 12.2. Restoring The Function Settings To The Factory Set Values

1. Turn the power OFF.
  2. Press and hold the **TARE** switch and press the **ZERO ON/OFF** switch to turn the power ON and to display the software version. Release the **ZERO ON/OFF** switch but continue to press the **TARE** switch until **Lrfno** with "no" blinking is displayed.
  3. Press the **SAMPLE** switch. **Lrfno** changes to **Lrfno** with "no" blinking.
  4. When **Lrfno** is displayed, press the **TARE** switch. The factory set values are restored. After displaying **End**, the scale goes to the weighing mode.
- To cancel the restoring procedure, press and hold the **ZERO ON/OFF** switch to turn off the scale.

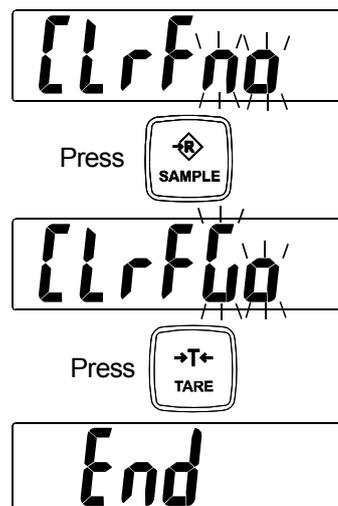

## 12.3.Function List

Item	Parameter	Description	
Key lock <i>Lock</i>	◆ 0	All function is enable.	
	1	Enable function: ON/OFF key, zero key, Tare key	
Auto power-off function <i>PoFF</i>	0	Auto power-off function disabled	Turns the power OFF automatically.
	◆ 1	Turns off after 5 minutes	
	2	Turns off after 10 minutes	
	3	Turns off after 15 minutes	
	4	Turns off after 30 minutes	
	5	Turns off after 60 minutes	
Weight display resolution <i>rESo</i>	0	1/3,000	Changes the minimum display.
	◆ 1	1/6,000 or 1/7,500	
	2	1/15,000 or 1/12,000	
	3	1/30,000	
Weighing unit <i>Unit</i>	SAMPLE	Proceeds to the next unit	Refer to "6. <b>SELECTING A WEIGHING UNIT</b> "
	TARE	Selects whether a unit is active or inactive	
	MODE	Proceeds to the next setting	
Zero tracking <i>trc</i>	0	Zero tracking function disabled	Tracks the zero drift.
	◆ 1	Zero tracking function enabled	
Weighing stability / response speed <i>Cond</i>	0		Response = Time from placing an object on the pan to turning on the stable indicator.
	1		
	◆ 2		
	3		
Backlight control <i>Lite</i>	0	Backlight always off	Sets the timing to turn off the backlight. Backlight turns on by weight change or switch operation.
	1	Backlight always on	
	◆ 2	Turns off 5 seconds after stabilizing	
	3	Turns off 10 seconds after stabilizing	
	4	Turns off 15 seconds after stabilizing	
	5	Turns off 30 seconds after stabilizing	
Brightness of backlight <i>L - i</i>	0	Dark	Adjusts the brightness of the backlight.
	1		
	◆ 2		
	3		
	4		
4	Bright		
Decimal point <i>Pnt</i>	◆ 0	Dot	
	1	Comma	

◆ Factory setting

Item	Parameter	Description	
Comparator mode <i>[P-L]</i>	0	Five-level (Result LED blinks.)	Sets comparator mode.
	1	Five-level (Result LED lights.)	
	2	Three-level (Result LED blinks.)	
	♦ 3	Three-level (Result LED lights.)	
	4	Seven-level (Result LED blinks.)	
	5	Seven-level (Result LED lights.)	
Comparison conditions <i>[P]</i>	0	Comparator disabled	Sets comparison conditions.  d = minimum display in kg.
	♦ 1	Compares all data	
	2	Compares all stable data	
	3	Compares all data of $\geq +5d$ or $\leq -5d$	
	4	Compares stable data of $\geq +5d$ or $\leq -5d$	
	5	Compares data of $\geq +5d$	
	6	Compares stable data of $\geq +5d$	
Comparator LED brightness <i>[P- ,]</i>	0	Dark	Adjusts LED brightness of comparison result.
	1	↑ ↓	
	♦ 2		
	3		
	4	Bright	
Normal/Negative comparison <i>[P-P]</i>	♦ 0	Normal comparison	Refer to "9. AUTO-TARE".
	1	Negative comparison for take-away check weighing	
Auto-tare function <i>At</i>	♦ 0	Auto-tare function disabled	
	1	Auto-tare function enabled	
Auto-tare timing <i>At-t</i>	0	Immediately after OK and stable	Timing to tare automatically after the comparison OK and stable weight.
	1	0.5 second after OK and stable	
	♦ 2	1.0 second after OK and stable	
	3	1.5 seconds after OK and stable	
	4	2.0 seconds after OK and stable	
	5	2.5 seconds after OK and stable	
	6	3.0 seconds after OK and stable	
	7	4.0 seconds after OK and stable	
	8	5.0 seconds after OK and stable	
	9	10 seconds after OK and stable	
Auto-tare of the initial weight <i>At-F</i>	♦ 0	Function disabled	Automatic operation.
	1	Tares the initial (container) weight	

♦ Factory setting

Item	Parameter	Description
Output mode <i>Print</i>	0	Stream mode / command mode
	1	Command mode only
	◆ 2	Output by print key / Command mode
	3	Auto print +/- data output / Command mode
	4	Auto print + data output / Command mode
	5	Auto comparator +/- data output on comparator OK / Command mode
	6	Auto comparator + data output on comparator OK / Command mode
		SJ-WP-BT model only


# 13. MAINTENANCE

## 13.1. Notes On Maintenance

- ❑ Do not disassemble the scale. Contact your local A&D dealer if the scale needs service or repair.
- ❑ Use the original packaging for transportation.
- ❑ Do not use organic solvents to clean the scale. Use a warm lint free cloth dampened with a mild detergent.
- ❑ Calibrate the scale periodically to maintain the weighing accuracy.

## 13.2. Pick cleaning

When unable to weigh properly due to dust between the load cell and overload stopper, insert a cleaning pick to remove the dust from the load cell.

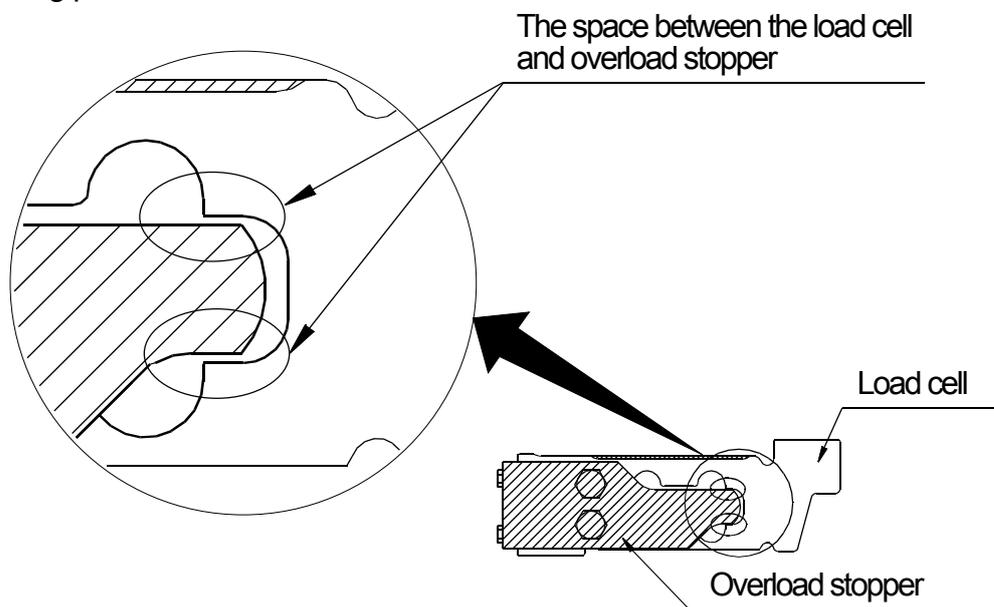


**Note:**

**There are five cleaning picks per sheet.**

**To use, cut off a cleaning pick when needed.**


**Use scissors or any bladed object to cut along the perforations to separate a cleaning pick for use.**

**Exercise extreme caution when using the bladed object to prevent personal injury.**

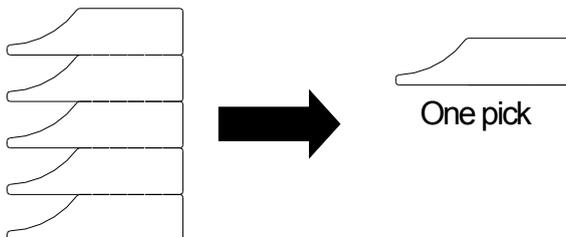

## 13.3. Error Codes

### Overload error

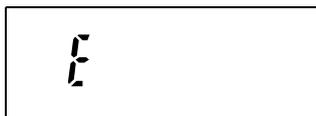

Indicates that an object beyond the weighing capacity has been placed on the weighing pan.


Remove the object from the weighing pan.

### Underload error



Indicates that the weight sensor receives a strong upward force. Check if there is anything sandwiched around the weighing pan. There is a possibility that the weight sensor or internal circuit may have a problem.

### Power-on zero error

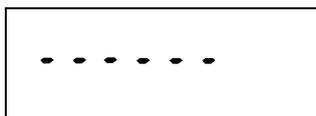


Indicates that the power is turned on with a load beyond the power-on zero range, or the weight value too unstable to perform power-on zero. Remove the load, or check if there is wind, vibration or anything touching the weighing pan.

### Unit weight error



Indicates that total weight of samples is too light to set the unit weight in the counting mode. Increase the number of samples and try again

### CAL error

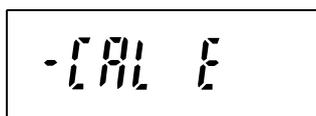

Indicates that the calibration procedure is canceled because the calibration weight is too light.

Check that the weighing pan is installed properly and the mass of the calibration weight is correct.


### Low battery



Indicates that the batteries have run out.

Replace them with new batteries.

### Other

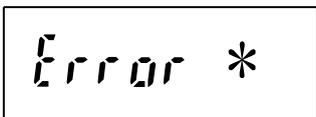

There may be an internal malfunction.

(\* indicates an error number.)

**Note:** If the error persists or other errors occur, contact your local A&D dealer.


# 14. SPECIFICATIONS

## Specifications

Model	SJ-3000WP/-BT	SJ-6000WP/-BT	SJ-15KWP/-BT	SJ-30KWP/-BT
Capacity	3 kg	6 kg	15 kg	30 kg
Minimum display "d"	0.001 kg	0.002 kg	0.005 kg	0.01 kg
	0.0005 kg*	0.001 kg*	0.002 kg*	0.005 kg*
	0.0002 kg	0.0005 kg	0.001 kg	0.002 kg
	0.0001 kg	0.0002 kg	0.0005 kg	0.001 kg
Counting	Number of samples: 5, 10, 20, 50 or 100 pieces Maximum count: 120,000			
Repeatability (SD)	0.5 g	1 g	2 g	5 g
Linearity	±1 g	±2 g	±5 g	±10 g
Sensitivity drift	±50 ppm / °C (5 °C to 35 °C / 41 °F to 95 °F)			
Display	Weight display: 7 segment LCD with Backlight, Character height: 26 mm Comparator LEDs: red / yellow / green / yellow / red			
Display update	20 times per second			
Operating temperature	-10 °C to 40 °C / 14 °F to 104 °F, Less than 85 %RH			
Power	6 x R20P / LR20 / "D" size batteries			
Battery life (Approximately) SJ-WP Model	5000 hours with alkaline cells at 20 °C (LED & Backlight off) 2000 hours with manganese cells at 20 °C (LED & Backlight off)			
Battery life (Approximately) SJ-WP-BT Model	2500 hours with alkaline cells at 20 °C (LED & Backlight off) 1300 hours with manganese cells at 20 °C (LED & Backlight off)			
Pan size	232 (W) x 192 (D) mm / 9.13 (W) x 7.56 (D) in.			
Dimensions	236 (W) x 260 (D) x 148 (H) mm / 9.3 (W) x 10.2 (D) x 5.8 (H) in.			
Mass	Approximately 4 kg / 9 lb			
Accessories	This manual, Screwdriver, Cleaning pick (One sheet)			
Sold separately	Cleaning pick (Five sheets) AXP-094038331			

\* : Factory setting

Unit: mm

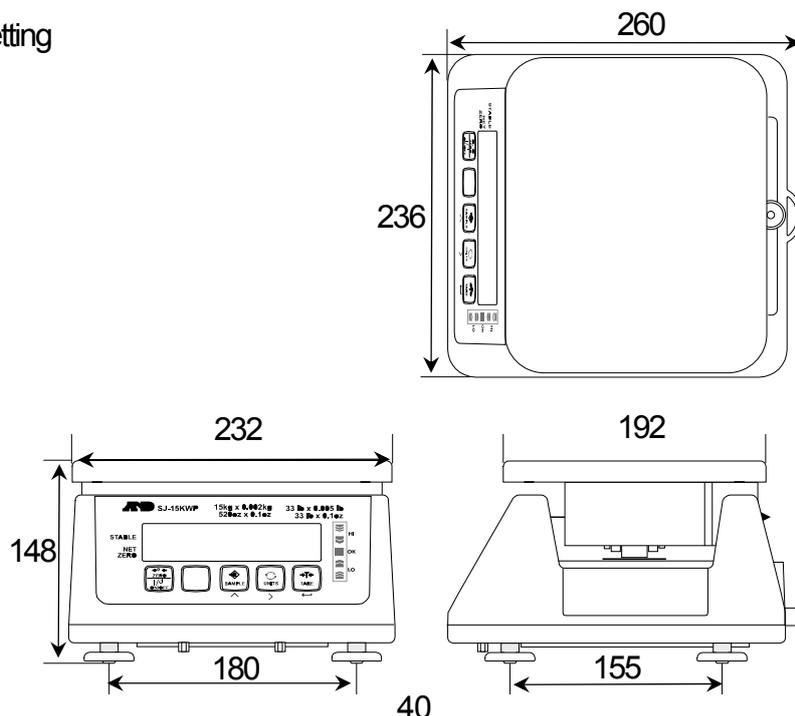

## Other weighing units

Model		<b>SJ-3000WP</b>	<b>SJ-6000WP</b>	<b>SJ-15KWP</b>	<b>SJ-30KWP</b>
g	Capacity	3000 g	6000 g	15000 g	30000 g
	Minimum display	1 g	2 g	5 g	10 g
		0.5 g*	1 g*	2 g*	5 g*
		0.2 g	0.5 g	1 g	2 g
		0.1 g	0.2 g	0.5 g	1 g
lb	Capacity	6.6 lb	13 lb	33 lb	66 lb
	Minimum display	0.002 lb	0.005 lb	0.01 lb	0.02 lb
		0.001 lb*	0.002 lb*	0.005 lb*	0.01 lb*
		0.0005 lb	0.001 lb	0.002 lb	0.005 lb
		0.0002 lb	0.0005 lb	0.001 lb	0.002 lb
oz	Capacity	105 oz	210 oz	520 oz	1050 oz
	Minimum display	0.05 oz	0.1 oz	0.2 oz	0.5 oz
		0.02 oz*	0.05 oz*	0.1 oz*	0.2 oz*
		0.01 oz	0.02 oz	0.05 oz	0.1 oz
		0.005 oz	0.01 oz	0.02 oz	0.05 oz
ozt	Capacity	96 ozt	193 ozt	480 ozt	960 ozt
	Minimum display	0.05 ozt	0.1 ozt	0.2 ozt	0.5 ozt
		0.02 ozt*	0.05 ozt*	0.1 ozt*	0.2 ozt*
		0.01 ozt	0.02 ozt	0.05 ozt	0.1 ozt
		0.005 ozt	0.01 ozt	0.02 ozt	0.05 ozt
lb-oz	Capacity	6 lb 9 oz	13 lb	33 lb	66 lb
	Minimum display	0.1 oz	0.1 oz	0.1 oz	0.1 oz
Catty-tl ( <i>HG</i> )**	Capacity	4 c 15 tl	9 c 14 tl	24c 12tl	49 c 9 tl
	Minimum display	0.01 tl	0.1 tl	0.1 tl	0.1 tl
Catty-tl ( <i>HJ</i> )**	Capacity	5 c	10 c	25 c	50 c
	Minimum display	0.01 tl	0.1 tl	0.1 tl	0.1 tl
Catty-tl ( <i>T</i> )**	Capacity	5 c	10 c	25 c	50 c
	Minimum display	0.01 tl	0.1 tl	0.1 tl	0.1 tl
Tola	Capacity	257 t	510 t	1280 t	2570 t
	Minimum display	0.1 t	0.2 t	0.5 t	1 t
		0.05 t*	0.1 t*	0.2 t*	0.5 t*
		0.02 t	0.05 t	0.1 t	0.2 t
		0.01 t	0.02 t	0.05 t	0.1 t

\* : Factory setting

\*\* : Catty-tael, *HG* : Hong Kong General / Singapore, *HJ* : Hong Kong Jewelry, *T* : Taiwan

## SJ-WP Series Class III Model Specifications

MODEL	SJ-3000WP/BT	SJ-6000WP/BT	SJ-15KWP/BT	SJ-30KWP/BT
Accuracy class	III			
Maximum Capacity	3000 g	6000 g	15 kg	30 kg
Readability	0.5 g(0~1500 g) 1 g(1500~3000 g)	1 g(0~3000 g) 2 g(3000~6000 g)	2 g(0~6 kg) 5 g(6~15 kg)	5 g(0~15 kg) 10 g(15~30 kg)
Minimum capacity	10 g	20 g	40 g	100 g
Maximum tare	3000 g	6000 g	15 kg	30 kg
Operating temp.	-10°C~40°C (No condensation)			
Display	Weight display: 7 segment liquid crystal display with White Backlight Character height: 26mm Comparator LEDs: red / yellow / green / yellow / red			
Display update	Approximately 20 times per second			
Power	6 x R20P / LR20 / "D" size batteries			
Battery life (Approximately) SJ-WP Model	5000 hours with alkaline cells at 20°C (LED & Backlight off ) 2000 hours with manganese cells at 20°C (LED & Backlight off )			
Battery life (Approximately) SJ-WP-BT Model	2500 hours with alkaline cells at 20°C (LED & Backlight off ) 1300 hours with manganese cells at 20°C (LED & Backlight off )			
Platform size	232 (W) x 192 (D) mm			
Dimensions	236 (W) x 260 (D) x 148 (H) mm			
Weight	Approximately 4 kg			
Calibration weight	3000 g ± 0.1 g	6000 g ± 0.2 g	15 kg ± 0.5 g	30 kg ± 1 g

### Note:

- The range for power-on zero is within ± 10% of the weighing capacity at the calibrated zero point.
- The weighing units "g", "kg" and "pcs" are available.
- The selection in the function settings rES0, AL, AL-t, and AL-F is not available.

### Power-on:

Turning the power on, the scale will be automatically set to zero.

If the power is switched on with a load within +/- 10 % of the weighing capacity at the calibrated zero point (power-on zero range), the scale is zeroed and the ZERO indicator turns on.

If the load is beyond the power-on zero range, the scale is tared and the ZERO and the NET indicators turn on.


### Error Message:

- E : Indicates that the gross value (weight value with no tare operation) is less than -19d. If the STABLE indicator is ON, press the ZERO  
ON/OFF switch to ZERO the scale. If the STABLE indicator is OFF, turn the power off and on again.
- If these instructions do not work, there is a possibility that the weight sensor or internal circuit may have a problem.


# 15. GRAVITY ACCELERATION

## Values of gravity at various locations

Amsterdam	9.813 m/s <sup>2</sup>	Manila	9.784 m/s <sup>2</sup>
Athens	9.807 m/s <sup>2</sup>	Melbourne	9.800 m/s <sup>2</sup>
Auckland NZ	9.799 m/s <sup>2</sup>	Mexico City	9.779 m/s <sup>2</sup>
Bangkok	9.783 m/s <sup>2</sup>	Milan	9.806 m/s <sup>2</sup>
Birmingham	9.813 m/s <sup>2</sup>	New York	9.802 m/s <sup>2</sup>
Brussels	9.811 m/s <sup>2</sup>	Oslo	9.819 m/s <sup>2</sup>
Buenos Aires	9.797 m/s <sup>2</sup>	Ottawa	9.806 m/s <sup>2</sup>
Calcutta	9.788 m/s <sup>2</sup>	Paris	9.809 m/s <sup>2</sup>
Cape Town	9.796 m/s <sup>2</sup>	Rio de Janeiro	9.788 m/s <sup>2</sup>
Chicago	9.803 m/s <sup>2</sup>	Rome	9.803 m/s <sup>2</sup>
Copenhagen	9.815 m/s <sup>2</sup>	San Francisco	9.800 m/s <sup>2</sup>
Cyprus	9.797 m/s <sup>2</sup>	Singapore	9.781 m/s <sup>2</sup>
Djakarta	9.781 m/s <sup>2</sup>	Stockholm	9.818 m/s <sup>2</sup>
Frankfurt	9.810 m/s <sup>2</sup>	Sydney	9.797 m/s <sup>2</sup>
Glasgow	9.816 m/s <sup>2</sup>	Taichung	9.789 m/s <sup>2</sup>
Havana	9.788 m/s <sup>2</sup>	Tainan	9.788 m/s <sup>2</sup>
Helsinki	9.819 m/s <sup>2</sup>	Taipei	9.790 m/s <sup>2</sup>
Kuwait	9.793 m/s <sup>2</sup>	Tokyo	9.798 m/s <sup>2</sup>
Lisbon	9.801 m/s <sup>2</sup>	Vancouver, BC	9.809 m/s <sup>2</sup>
London (Greenwich)	9.812 m/s <sup>2</sup>	Washington DC	9.801 m/s <sup>2</sup>
Los Angeles	9.796 m/s <sup>2</sup>	Wellington NZ	9.803 m/s <sup>2</sup>
Madrid	9.800 m/s <sup>2</sup>	Zurich	9.807 m/s <sup>2</sup>

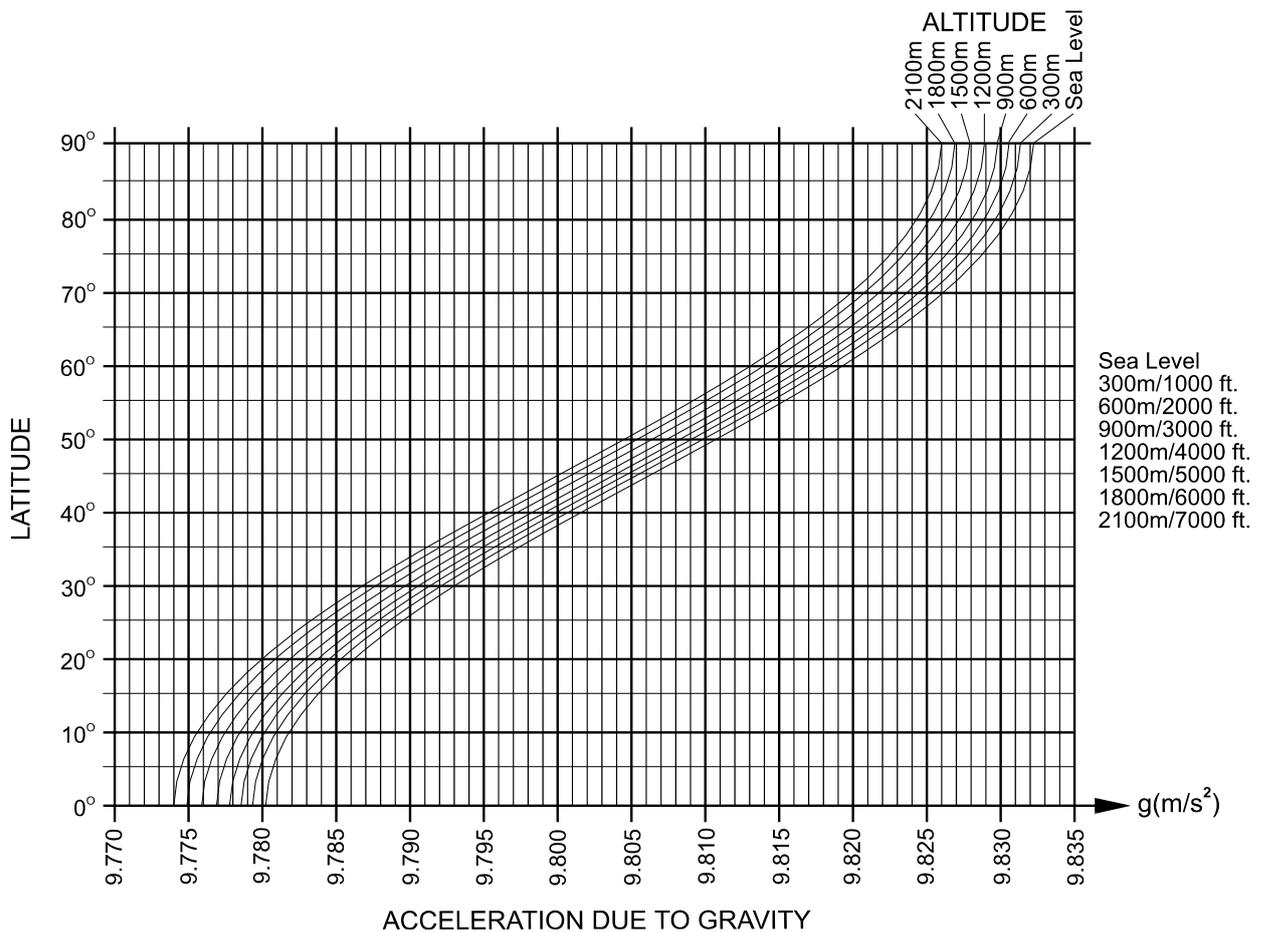

World map

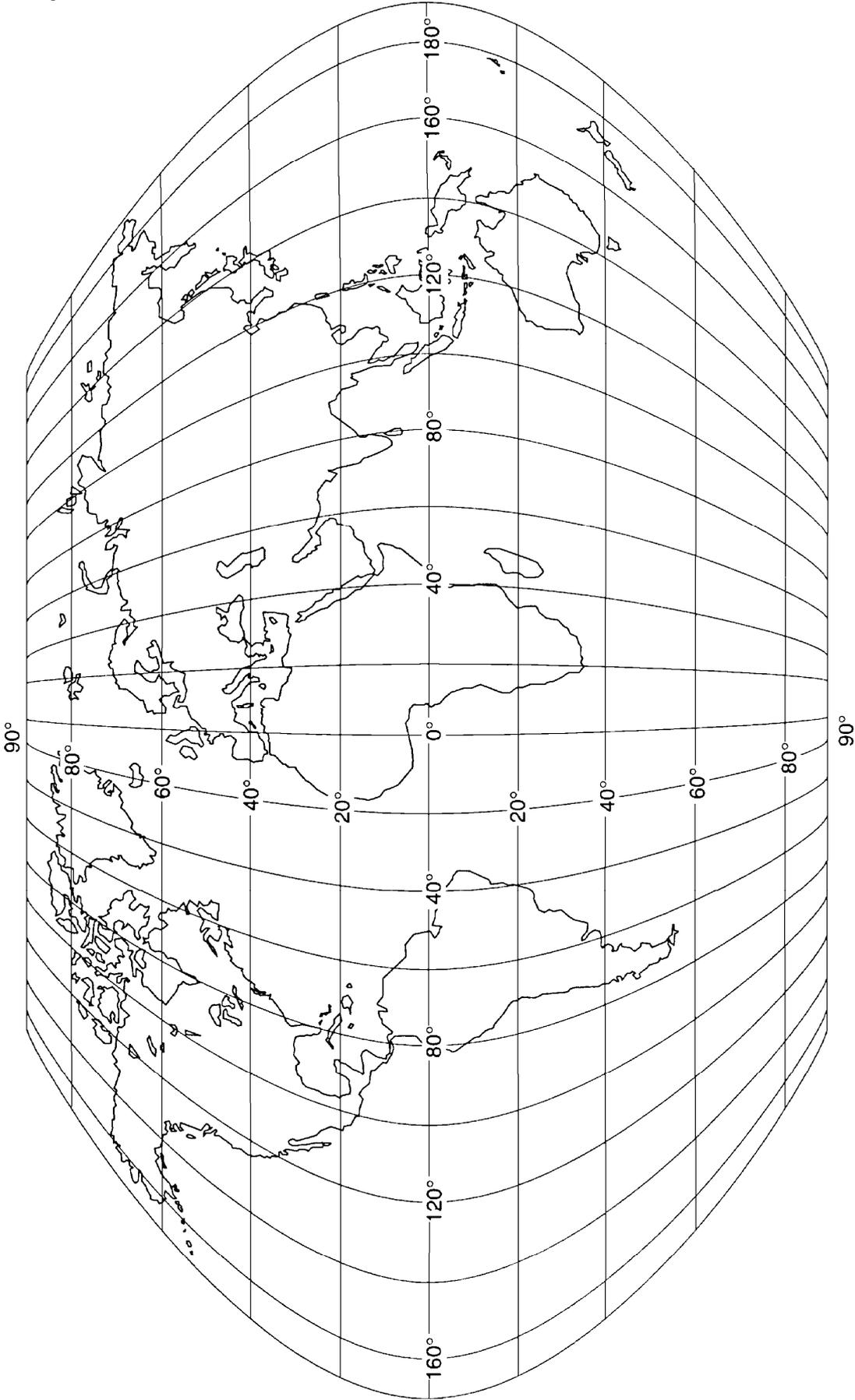





#### **A&D Company, Limited**

3-23-14 Higashi-Ikebukuro, Toshima-ku, Tokyo 170-0013, JAPAN  
Telephone: [81] (3) 5391-6132 Fax: [81] (3) 5391-1566

#### **A&D ENGINEERING, INC.**

1756 Automation Parkway, San Jose, California 95131, U.S.A.  
Telephone: [1] (408) 263-5333 Fax: [1] (408) 263-0119

#### **A&D INSTRUMENTS LIMITED**

Unit 24/26 Blacklands Way, Abingdon Business Park, Abingdon, Oxfordshire OX14 1DY United Kingdom  
Telephone: [44] (1235) 550420 Fax: [44] (1235) 550485

#### **A&D AUSTRALASIA PTY LTD**

32 Dew Street, Thebarton, South Australia 5031, AUSTRALIA  
Telephone: [61] (8) 8301-8100 Fax: [61] (8) 8352-7409

#### **A&D KOREA Limited**

한국에이.엔.디(주)  
서울특별시 영등포구 국제금융로6길33 (여의도동) 맨하탄빌딩 817 우편 번호 07331  
( 817, Manhattan Bldg., 33. Gukjegeumyung-ro 6-gil, Yeongdeungpo-gu, Seoul, 07331 Korea )  
전화: [82] (2) 780-4101 팩스: [82] (2) 782-4264

#### **ООО A&D RUS**

ООО "ЭЙ энд ДИ РУС"  
121357, Российская Федерация, г.Москва, ул. Верейская, дом 17  
( Business-Center "Vereyskaya Plaza-2" 121357, Russian Federation, Moscow, Vereyskaya Street 17 )  
тел.: [7] (495) 937-33-44 факс: [7] (495) 937-55-66

#### **A&D INSTRUMENTS INDIA PRIVATE LIMITED**

ऐ&डी इन्स्ट्रूमेन्ट्स इण्डिया प्रा० लिमिटेड  
509, उद्योग विहार, फेस -5, गुडगांव - 122016, हरियाणा, भारत  
( 509, Udyog Vihar, Phase-V, Gurgaon - 122 016, Haryana, India )  
फोन : 91-124-4715555 फैक्स : 91-124-4715599