

Professional Blood Pressure Monitor

Model: UM-301

General Specifications

- **Energy-efficient**

The mechanical exhaust valve significantly reduces the number of battery replacements required. By incorporating a mechanical exhaust valve into the cuff tube, it can be used continuously for approximately 120 days when performing 50 measurements per day. One battery can perform approximately 6,000 measurements.

- **Interchangeable cuffs**

Mechanical design allows for interchangeable cuff sizes.

- **Slim & portable**

- **Flexible grip**

A robust yet flexible grip accommodating diverse gripping styles.

- **One-touch use**

Intuitive and simple single-tap operation to prevent mistakes.



SWANHEART Grip

*"Easy to Use.
Built to Last."*



*The reliable one-hand
blood pressure monitor for
medical professionals.*



Mechanical Exhaust Valve

Products with built-in electric exhaust valves use up batteries quickly due to electronic control.

The UM-301 has a built-in mechanical exhaust valve in the cuff tube. This mechanical exhaust valve does not use electricity, so it can be used for long periods of time without battery replacement.

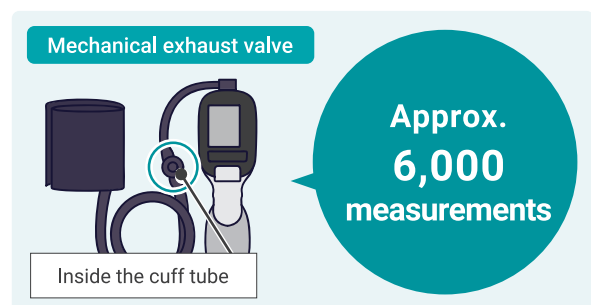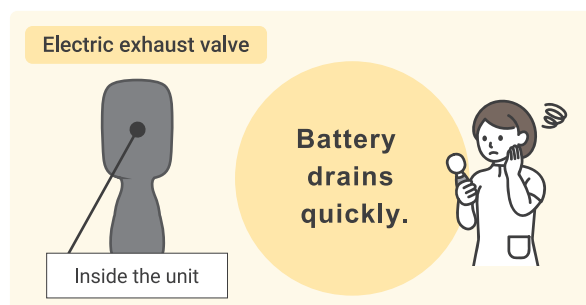


Interchangeable Cuffs

Mechanical design allows for interchangeable cuff sizes. Mechanical exhaust valves matched to each cuff size are included with the cuffs.



For a more comfortable daily routine.

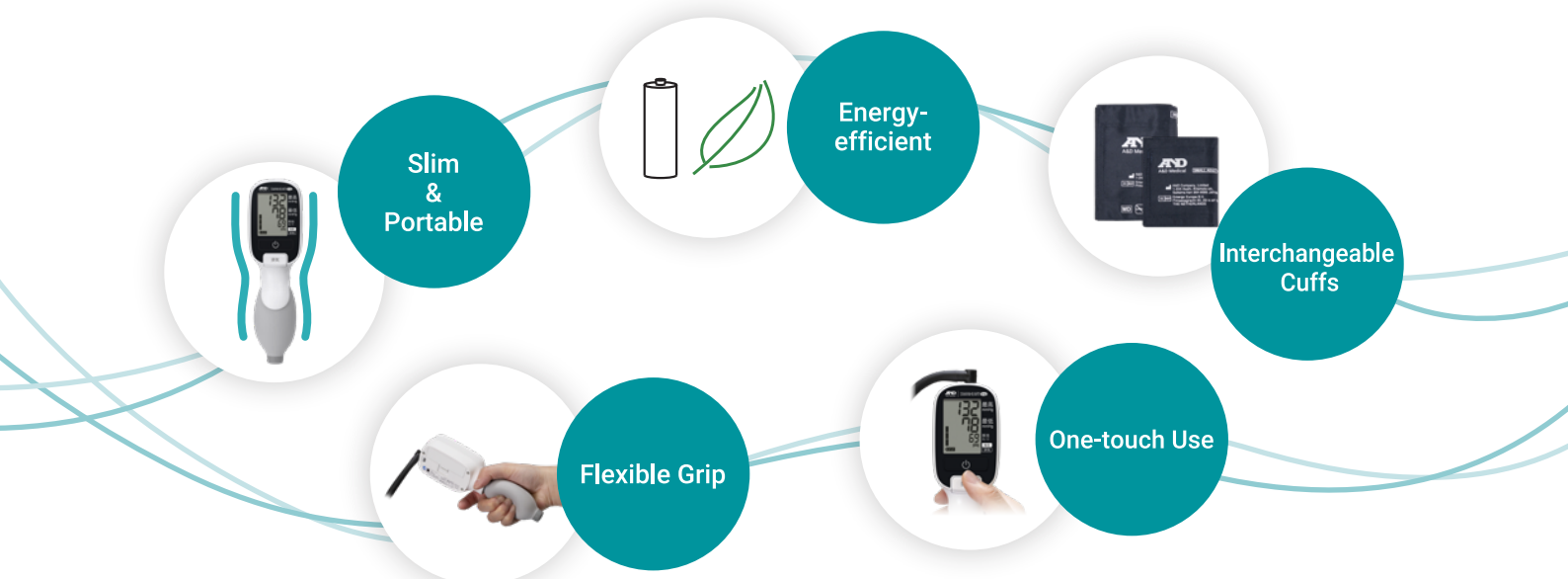

Technical Specifications

Measurement Method	Oscillometric
Measurement Range	Pressure: 0 to 299 mmHg / Pulse: 40 bpm to 180 bpm SYS: 60 mmHg to 279 mmHg / DIA: 40 mmHg to 200 mmHg
Measurement Accuracy	Pressure: ± 3 mmHg / Pulse: ± 5 % of the reading value
Power Source	AA alkaline (or manganese) battery x 1
Battery Life	Approx. 6,000 measurements with alkaline battery / Approx. 2,000 measurements with manganese battery <small>*Based on continuous measurements at room temperature (23 °C) with each measurement taking 50 seconds.</small>
Dimensions (W x H x D)	Approx. 53 x 42 x 182 mm
Unit Weight	Approx. 120 g excl. battery

Optional Accessories

Product Code	Option Name	Arm Circumference
UM-AURS6K6KJC	A Size Cuff	22 cm - 32 cm
UM-SARS6K6KJC	SA Size Cuff	16 cm - 24 cm
AX-PA4003048	Inflation Bulb	-



UM-AURS6K6KJC



UM-SARS6K6KJC



AX-PA4003048

What Is In The Box

- 1 Blood Pressure Monitor
- 1 A Size Cuff (22 cm - 32 cm)
- 1 AA Battery
- 1 Instructions For Use

Specifications may differ in each country.
Please contact A&D subsidiary for more information.

A&D
A&D Medical

A&D Company Limited (JAPAN)
URL: aandd.jp

A&D ENGINEERING, INC. (USA)
URL: andonline.com

A&D Instruments Canada Inc. (CANADA)
URL: autocontrol.com

A&D INSTRUMENTS LIMITED (United Kingdom)
URL: andprecision.com

A&D AUSTRALASIA PTY LTD (Australia)
URL: andaustrolasia.com.au

A&D RUS CO., LTD. (Russia)
URL: and-rus.ru

A&D Technology Trading (Shanghai) Co., Ltd (China)
URL: aanddtech.cn

A&D INSTRUMENTS INDIA PRIVATE LIMITED (India)
URL: aanddindia.in

A&D Instruments (Thailand) Limited
URL: thai.andprecision.com/

A&D KOREA Limited
URL: andk.co.kr/main/