

# Accessories for Micro-Analytical Balances

**AND**  
A&D Weighing

## For cleaning balances:

### AX-CLEANING-SET

Balance Cleaning Set

- Includes two sizes of bristle brushes and a chamois leather cloth.
- The bristles are made of animal fiber, which has excellent antistatic performance.
- The soft chamois leather is gentle on the balance surface while efficient at trapping dirt.
- Highly convenient and recommended for daily maintenance of balances.

Material (Handle): Wood  
Material (Bristles): Goat Hair  
Brush Length: Approx. 220 mm  
Brush Width: (Large Brush) 27 mm  
(Small Brush) 14 mm



## For handling calibration weights:

### AX-TWEEZERS-25

Tweezers, Long type, 25 cm



- Tip shape designed to handle various sizes of calibration weights (i.e. plate weights of 1 to 500 mg, small weights of 1 to 5 g, and large weights of 10 to 500 g)
- Long enough to reach the weighing pan without having to put a hand inside the weighing chamber and thereby causing temperature disturbance.

### AX-TWEEZERS-13B

Tweezers, Bent-tip type, 13 cm



### AX-TWEEZERS-14S

Tweezers, Straight-tip type, 14 cm



- Shorter length suitable for micro-analytical balances with a small weighing chamber, such as the BA-T/BA series.

## For holding test tubes, measuring flasks, etc.:

### AX-HOLDER-SET-B

Sample Cup Holder Set (Large/Medium/Small)

- The holders can be placed on a balance and used to secure a tilted container for measurement.
- The holders can also be turned over and used to secure a long test tube.

Material: Conductive antistatic silicone  
Large size : 67 (W) × 67 (D) × 40 (H) mm, 45 g  
Medium size: 56 (W) × 56 (D) × 36 (H) mm, 32 g  
Small size: 40 (W) × 40 (D) × 32 (H) mm, 16 g



### BA-032

Microtube Holder

- Secures a microtube inside the BA-T/BA microbalance.
- Easy to adjust the angle for ergonomic use.

Applicable models: BA-6/BA-6D/BA-6TE/BA-6DTE  
Compatible microtubes: Outer diameter of 11 mm or less



## For weighing small amounts of samples:

### AX-TRAY-10ML

Disposable Weigh Boat

- Anti-static treatment to prevent static electricity.

Capacity: 10 ml  
Weight: 0.49 g  
Quantity: 250 pieces




\* Not suitable for use with microbalances

### AX-SPATULA

Micro Spatula

- Suited to handling especially tiny amounts of fine powders.

Length: 182 mm  
Weight: 4.3 g  
Quantity: 5 pieces

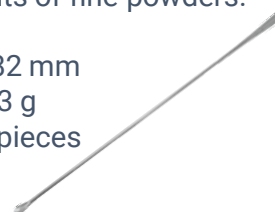

## Discover Precision

A&D Company, Ltd. (JAPAN)  
URL: aandd.jp  
A&D Engineering, Inc. (USA)  
URL: andonline.com  
A&D Australasia Pty Ltd. (Australia)  
URL: andaustalasia.com.au  
A&D Instruments Ltd. (United Kingdom)  
URL: andprecision.com  
<German Sales Office>  
URL: andprecision.com

A&D Korea Ltd. (South Korea)  
URL: andk.co.kr  
A&D Rus Co., Ltd. (Russia)  
URL: and-rus.ru  
A&D Instruments India (P) Ltd. (India)  
URL: aanddindia.in  
A&D Sciencetech Taiwan Ltd. (Taiwan)  
URL: aandd.com.tw  
A&D Instruments Thailand Ltd. (Thailand)  
URL: thai.andprecision.com